

Please keep this manual for future reference.



# SHIATSU FOOT & CALF MASSAGER

with Heat & RockerTech™

## User Manual



Infinity Massage Chairs  
72 Stard Road  
Seabrook, NH 03874  
603-910-5000

[www.infinitymassagechairs.com](http://www.infinitymassagechairs.com)

Thank you for choosing the Infinity Shiatsu Foot & Calf Massager with Heat & RockerTech™. Please read the User Manual carefully before using the massager. Keep this booklet in a safe place for future reference.


### **IMPORTANT SAFETY WARNINGS:**

- Please use this product in accordance with the guidelines of this manual.
- Do not allow any foreign objects to enter the clearance between the massage head and casing.
- To prevent short-circuiting, avoid getting the massager wet. Avoid contact with liquids/ water.
- If the main cable or plug is damaged, do not turn on the massager or attempt to repair the plug. Please contact Infinity for troubleshooting instructions.
- Disconnect the product from the power supply after use and before cleaning to avoid electrical injury.
- Children and disabled individuals should be supervised while using this product.
- Do not attach anything else to the massager during use.
- Do not place heavy objects on massager when not in use.
- Do not use this product within one hour after eating.
- Do not use this product immediately after there has been a sharp change in the environmental temperature.
- Do not use this product in a dusty or corrosive environment.
- Do not cover the massager with a blanket or pillow during use.
- Do not fall asleep while using this product.
- Do not use product on pets.

### **The following people should consult their doctor prior to using this massager:**

- People with contagious skin conditions.
- People with osteoporosis or other bone and muscle conditions.
- People with heart disease and/or people with electrical medical instruments such as pacemakers.
- People with high fevers.
- Women who are pregnant.
- People with injuries or recovering from injuries.
- The massager has a heated surface. Those sensitive to heat should use this massager with caution.

## MASSAGER OPERATION

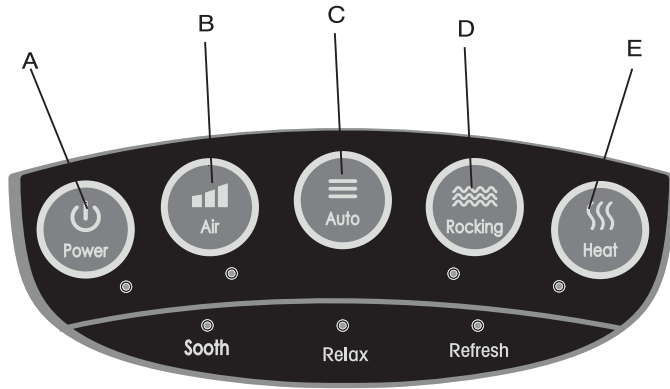

### A. Power button: Turns massager on/off

- Indicator light stays on when the massager is on. Power On defaults to the first Auto program.
- Once the massager is on, press the power button again to stop all functions. Indicator light will blink.

### B. Air button: Controls air pressure intensity

- Adjust the intensity level or turn airbags off with this button.
  - Mode A** (Indicator light blinks slowly): Low intensity
  - Mode B** (Indicator light blinks quickly): Medium intensity
  - Mode C** (Indicator light stays on): High intensity
  - Mode D** (No light): Air off
- Once a mode is selected, press the intensity button again to change the level of intensity.

### C. Auto button: Controls Auto programs

- Cycle through the three available Auto programs with this button.
  - 1) Soothe:** Rolling message nodes are inactive while airbags and RockerTech™ are activated.
  - 2) Relax:** Airbags, rolling message nodes, and RockerTech are activated. Rolling nodes work intermittently.
  - 3) Refresh:** Airbags, rolling message nodes, and RockerTech are activated. Rolling nodes work continuously.
- Press Auto a fourth time to enter Manual mode, in which you can manually adjust Air, RockerTech, and Heat settings.

### D. Rocking button: Controls RockerTech™ tapping rhythm

- Select from three RockerTech modes with this button.
  - Mode A** (Indicator light stays on): RockerTech is continuous.
  - Mode B** (Indicator light blinks quickly): RockerTech functions in varied pattern.
  - Mode C** (Indicator light blinks slowly): RockerTech is intermittent.
- Once an auto program is selected, press the Rocking button again to stop the current RockerTech mode. To select another one, simply cycle through with the Rocking button.

### E. Heating button: Turns on heating function

- Indicator light will turn on when foot heating function is turned on. Light will turn off when heat is off.
- Once heating function is on, press the heating button again to deactivate the heat.
- Please use the heat function with caution. Do not keep heat on for prolonged periods of time.

### Massage Session Time & Auto Shut-off

- The massager will automatically end the massage session after 15 minutes. If another session is desired, please make a new selection.
- Multiple back-to-back uses could result in auto shut-off to prevent overheating. If this occurs, please wait at least 30 minutes to allow the motors to cool down before attempting to start a new session.

## MAINTENANCE

### A. Storage

- When not using the massager, store it back in the box and place it in a dry, well-ventilated, cool area.
- Do not twist wires.
- Avoid sharp objects to prevent scratching surface of the product.
- Keep out of direct sunlight and high temperatures.

### B. Cleaning

- Unplug power cord before cleaning.
- Use a dry cloth only to clean the surface of the product. Avoid direct contact with water.
- Do not use corrosive cleaners such as benzene or bleach products.
- The fabric pouch can be removed for cleaning.

### C. Maintenance

- Keep the product clean and away from moisture.
- Cover up the product when not in use to avoid dust.



**Warning:** This product has an internal temperature protection device. When the safe temperature or safe heat use duration is exceeded, the heat will be turned off automatically. Please unplug the power to allow machine to cool after use.

### Infinity Technical Support

If you have any questions about maintenance or are experience technical difficulties with your massager, please contact Infinity's technical support team at 603-910-5000 for assistance.

## TROUBLESHOOTING

Problem	Possible Cause	Possible Solution
Massager not responding	No power	Check power supply
	Function keys are not responding	Restart massager using On/Off button or unplug the device
Massager suddenly stops operating	No power	Check power supply
	Massager in-use for more than 15 minutes	Turn off massager to allow for cool-down period. Restart it using On/Off button.
	Overheating due to overuse	Turn off massager to allow for cool-down period. Restart it using On/Off button.
	Auto shut-off due to excessive force	Relax feet when massager is on.

## SPECIFICATIONS

<b>Model</b>	SHIATSU FOOT & CALF MASSAGER with Heat & RockerTech™
<b>Dimensions</b>	15.75" x 14.6" x 18.5" (L x W x H)
<b>Weight</b>	16.3 lbs
<b>Rated Voltage</b>	DC12V
<b>Rated Power</b>	60W
<b>Rated Duration</b>	15 Minutes