



Overture Massage Chair

User's Manual



Infinity Massage Chairs
72 Stard Road
Seabrook, NH 03874
603-910-5000
www.infinitymassagechairs.com

Please read SAFETY INSTRUCTIONS carefully before using the chair.

INFINITY

Overture Massage Chair

Function Introduction

1. Equipped with a patented 4D flexible swing massage mechanism on the backrest, the Overture massage chair is designed with functions such as kneading, tapping, shiatsu, knocking, and quick rubbing. This unique massage mechanism is able to creatively imitate human hand kneading techniques for an incredibly soothing massage.
2. Designed with heated ceramic massage rollers, the Overture provides heat therapy on the neck, shoulders, waist, and buttocks, relieving pain and promoting blood circulation.
3. Thanks to the chair's super long L-shaped track, the rollers are able to move from the neck area down to the thighs, providing a true full-body massage.
4. Zero wall space-saving technology allows the chair to be placed right up against the wall, as it will slide forward on its track during the massage and move back when finished. This means the Overture can be used in more compact spaces.
5. The zero gravity function can effectively relieve pressure on the spine and joints. This position places your legs above your heart, providing more oxygen to the brain and relaxing you both physically and mentally.
6. The Overture's convenient wireless full-color LCD screen remote control can display your real-time massage status and massage techniques. You can also use the remote to personalize your own custom massage program.
7. Bluetooth® control technology allows you to connect your smartphone, tablet, or other Bluetooth-enabled device and easily control your massage wirelessly through the Apple® or Android™ mobile app.
8. You can also enjoy your favorite music, audio books, or ambient soundscapes using your device, thanks to the Bluetooth-connected speakers in the headrest.
9. Restore strength and elasticity and reduce pain in your lower legs with the rolling calf massagers.
10. Relieve sore, tired feet with comprehensive foot units featuring rejuvenating foot rollers, gentle infrared heat, and airbag massage.
11. Airbag massage for the arms and hands help improve blood circulation and relaxation.
12. Multiple sets of airbags are also located in the waist and hip regions, providing total coverage to the waist, hip, and thighs with alternating air pressure techniques to loosen muscles.
13. The auto-extending footrest can be extended up to 8.5 inches.
14. Feel refreshed and energized with the Overture's air ionizer, which absorbs harmful substances in the air and improve air quality.
15. Easily control your massage with intelligent voice activation technology for a totally hands-free massage!



- ★ Thank you for choosing the Infinity Overture Massage Chair! Please read this manual carefully before using your new massage chair, and keep it in a safe place for future reference. Please also pay attention to the Safety Instructions.
- ★ Find us on YouTube at www.youtube.com/infinitymassagechairs for helpful tutorials and additional information on assembling and using your massage chair!
- ★ Please Note: We reserve the right to design modification. It is subject to change without notice. The pictures in this manual are only for reference. Please refer to the actual product, or contact Infinity Massage Chairs at 603-910-5000 for assistance.

Table of Contents





Preparation	Safety instructions	1
	External structure parts	5
	Internal structure parts	6
	Mechanical components diagram	7
	Mechanical components list	8
Installation	Installation	9
Massage Operation	Setting up your massage chair	15
	Remote Control	17
	Operation instructions & Function Illustration	23
	After massage	26
Maintenance	Cleaning and maintenance	27
Others	Malfunction guide	28
	Technical data	29

Please keep this manual in a safe place for future reference.

- ✦ Please read this manual carefully to ensure correct installation and operation of your massage chair.
- ✦ Safety Instructions are indicated as either Warning or Caution.

	Warning	This symbol indicates warning of serious danger.
	Caution	This symbol indicates caution of possible injury or damage to people or objects.




Symbol meaning






	Prohibited behavior.		Do not disassemble or attempt repair
	Must follow instructions to operate		Correct operation of cable




- ✦ Please keep the following information for future reference.




Warning

	<ul style="list-style-type: none"> • Those with osteoporosis, spinal injuries, or back problems • Those with a medical instrument embedded in the body, such as a pacemaker • Pregnant women • Those with skin disease or skin injuries • Those with a high fever • Intoxicated people • Cancer patients, heart disease patients, and acute care patients 	People with any of the indicated conditions should consult a doctor prior to using a massage chair.
	<ul style="list-style-type: none"> • Make sure there are no kids, pets, or other objects under, behind, or in front of the chair during operation. • If you feel any pain or discomfort during massage chair operation, stop use immediately and consult your physician. • Avoid massage on the stomach or knees, and use caution when massaging the neck. • Children and pets should not stand on the seat, back, or armrests of the chair. • Children and people with mental disabilities should only use the massage chair under supervision. 	Accident/injury
	<ul style="list-style-type: none"> • Insert the chair's plug into a suitable power outlet with the correct voltage. Ensure the plug is completely inserted into the socket to reduce risk of short circuit or fire. • Unplug the power cord when not in use. Stop all massage functions and turn off ON/OFF key before unplugging. • Do not operate the controller or massage chair with wet or damp hands. 	Fire and electric shock
	<ul style="list-style-type: none"> • Ensure there are no objects or debris blocking ventilation on any part of the massage chair. • The massage chair cushions are not meant to be folded and stored. 	Accident
	<ul style="list-style-type: none"> • Do not use the chair if you notice any damage (i.e.: split, knotted, twisted) to the power cord or plug. Do not plug into a loose or faulty outlet. • Only use the rated voltage recommended to avoid fire, electric shock, and damage. 	Fire and electric shock
	<ul style="list-style-type: none"> • If you experience any issues or dysfunctional behavior with the chair, or if it undergoes water damage, please contact Infinity Massage Chairs for repair and maintenance by a professional service agent. • If the power cord or plug is damaged, please contact Infinity directly for replacement. Please do not replace the power cord yourself. • In the case of any damage to the leather or cloth material of the chair, please stop using and contact Infinity for repair. 	Disassembly and maintenance


 Caution		
	<ul style="list-style-type: none"> • Please operate the chair on a flat or level surface. 	To avoid sudden falling.
	<ul style="list-style-type: none"> • Position the chair 4 inches from the wall to allow for sufficient space to recline up and down. 	To avoid damage to the chair or nearby walls
	<ul style="list-style-type: none"> • The maximum suggested operating time is 20 minutes. A long, continuous massage in the same position may cause some pain or discomfort. 	To prevent pain or injury from extensive massage time
	<ul style="list-style-type: none"> • Before using the chair, be sure to inspect everything to confirm that the massage rollers are in the correct position. 	To avoid accidental injury.
	<ul style="list-style-type: none"> • Be sure to store the remote control in the pocket when not in use to avoid dropping or damaging the remote. 	
	<ul style="list-style-type: none"> • Make sure that the power switch is in the OFF position before plugging the chair into the outlet. 	
	<ul style="list-style-type: none"> • Turn off all controls before removing the plug from the outlet. 	To avoid accident or injury when powering on
	<ul style="list-style-type: none"> • Make sure the plug is completely inserted into the socket.  	To reduce short circuits and any fire hazard.
	<ul style="list-style-type: none"> • Do not press or step on the massage rollers. • Do not place hands between the working massage heads. • Do not put hands or feet between mechanical parts. • Do not drop or insert anything into any openings or gaps of the massage chair. • Do not stand, sit, or place anything on the backrest, armrest, or legrest. • Do not stand on the chair during operation. • Do not move the chair during operation. • Do not pull out the plug or turn off the power during operation • Do not put your hands or head between or beneath the legrest. • Do not use the massage chair while wearing hair ornaments. • Do not keep hard items in pockets during thigh and buttocks massage. • Don't sit or press on the backrest when the chair is reclining/ reclined. 	To avoid accidents or injury
	<ul style="list-style-type: none"> • Keep the chair away from extreme heat, moisture, sharp objects, and corrosive or flammable products. • Do not use electric blankets or other heating products in the massage chair. • Do not place the massage chair in an environment with extreme amounts of dust or gas. • Do not place the massage chair in very small spaces or use in areas with poor air ventilation. 	To avoid fading, discoloring, or deterioration on upholstery or PU
	<ul style="list-style-type: none"> • Keep the chair away from stove, oven, or other heating appliances. • Do not place the chair directly under the sun. • Do not use the chair outdoors. 	

 Caution		
	<ul style="list-style-type: none"> • Don't use or store the chair in a wet environment, such as a pool or bathroom. • Keep the chair and remote control dry at all times. 	To avoid leakage or electric shock
	<ul style="list-style-type: none"> • Unplug the power before cleaning the chair. • Unplug the power after using the chair. 	To avoid electrical shock or fire hazard
	<ul style="list-style-type: none"> • Unplug the power if you see any damage on the plug or cord. 	
	<ul style="list-style-type: none"> • Hold the plug to unplug the power. Do not pull the wire to unplug the power. 	To avoid accidents
	<ul style="list-style-type: none"> • Please turn off the power switch and unplug the power after your massage. 	
	<ul style="list-style-type: none"> • In the event of a loss of electricity, unplug the massage chair's power cord immediately. 	To avoid damage or accidents when powering the chair on again


External structure

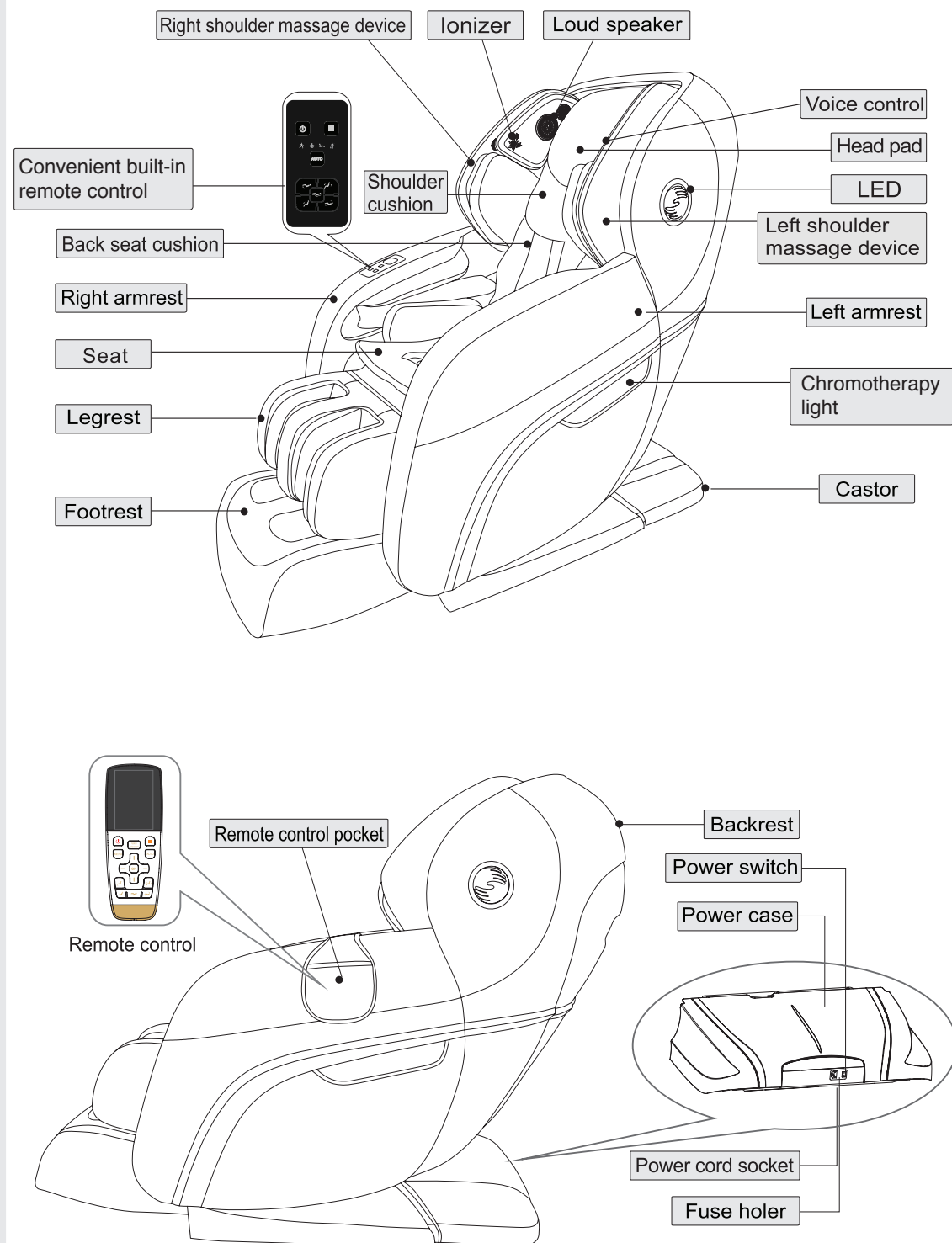

Internal structure

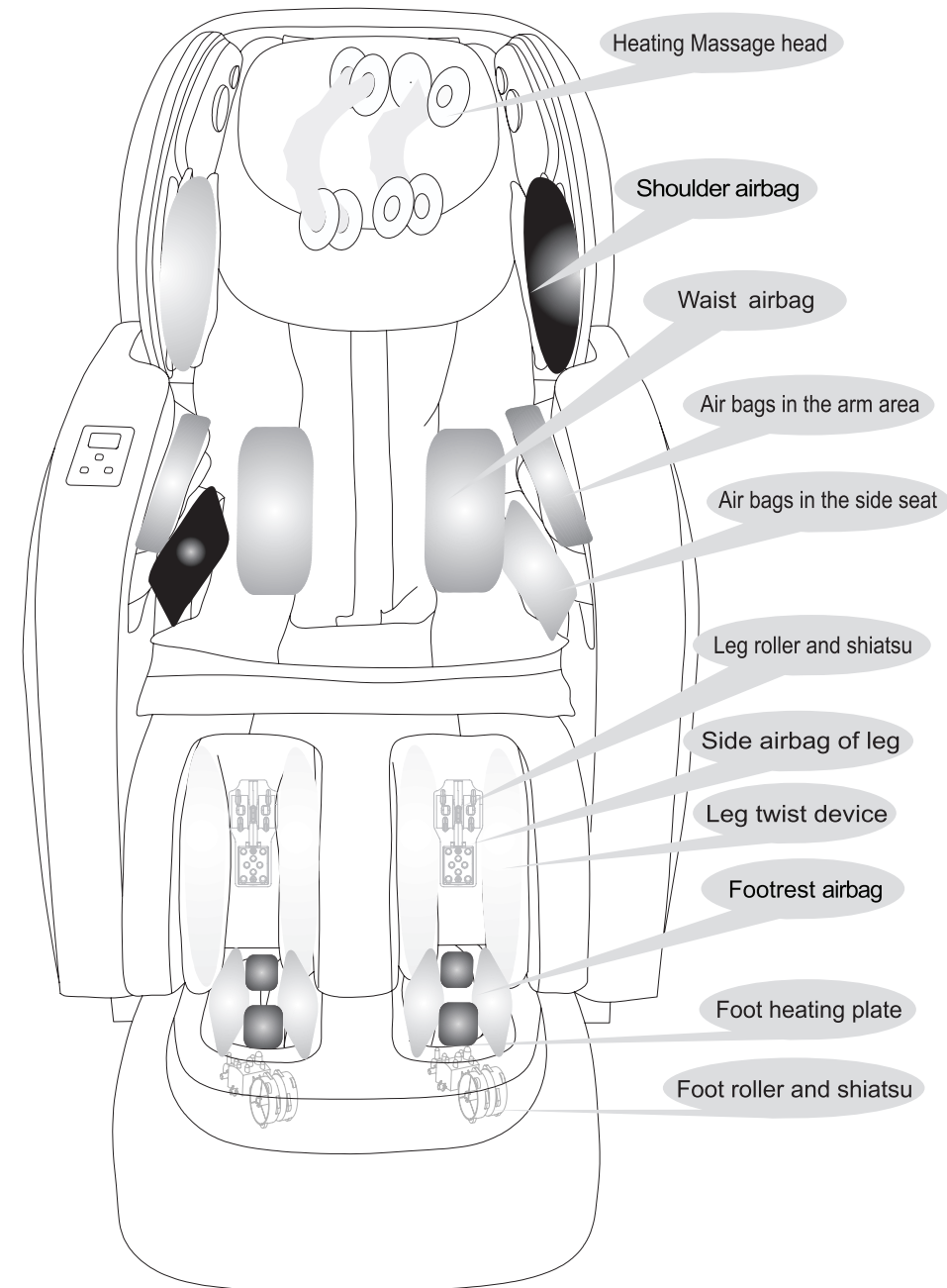

Explosion Diagram

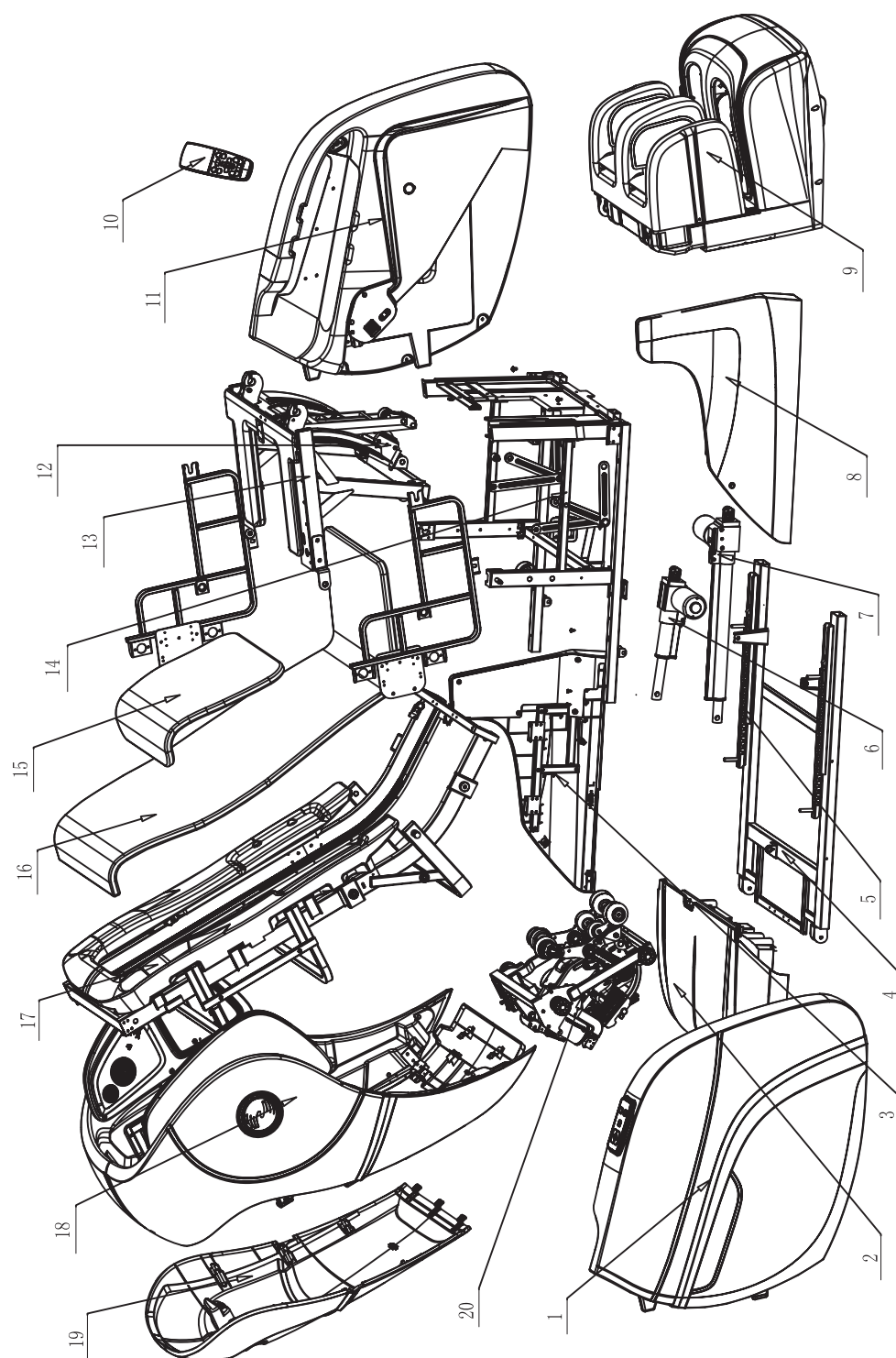

Preparation

Mechanical components list

NO.	Name	Q' ty	NO.	Name	Q' ty
1	Right armrest assembly	1			
2	Power case assembly	1			
3	Left side fender	1			
4	Base assembly	1			
5	Sliding rail	2			
6	Seat electric actuator	1			
7	Backrest electric actuator	1			
8	Right side fender	1			
9	Foot and leg assembly	1			
10	Remote control	1			
11	Left amrest assembly	1			
12	Footrest electric actuator	1			
13	Seat frame surface assembly	1			
14	Seat bottom frame assembly	1			
15	Pillow	1			
16	Back and seat cushion	1			
17	Backrest assembly	1			
18	Shoulder device assembly	1			
19	Back cover assembly	1			
20	Massage machine	1			


Preparation

Installation (1)

Step 1: Unpack the massage chair body and accessories from the box.

1. Remove the accessories boxes, chair body, backrest pad, pillow, power cord, remote control, and remote control pocket from the box.
2. Confirm all the necessary accessories are present.

Please Note: You can also visit our YouTube channel at www.youtube.com/infinitymassagechairs to watch our easy assembly tutorial!

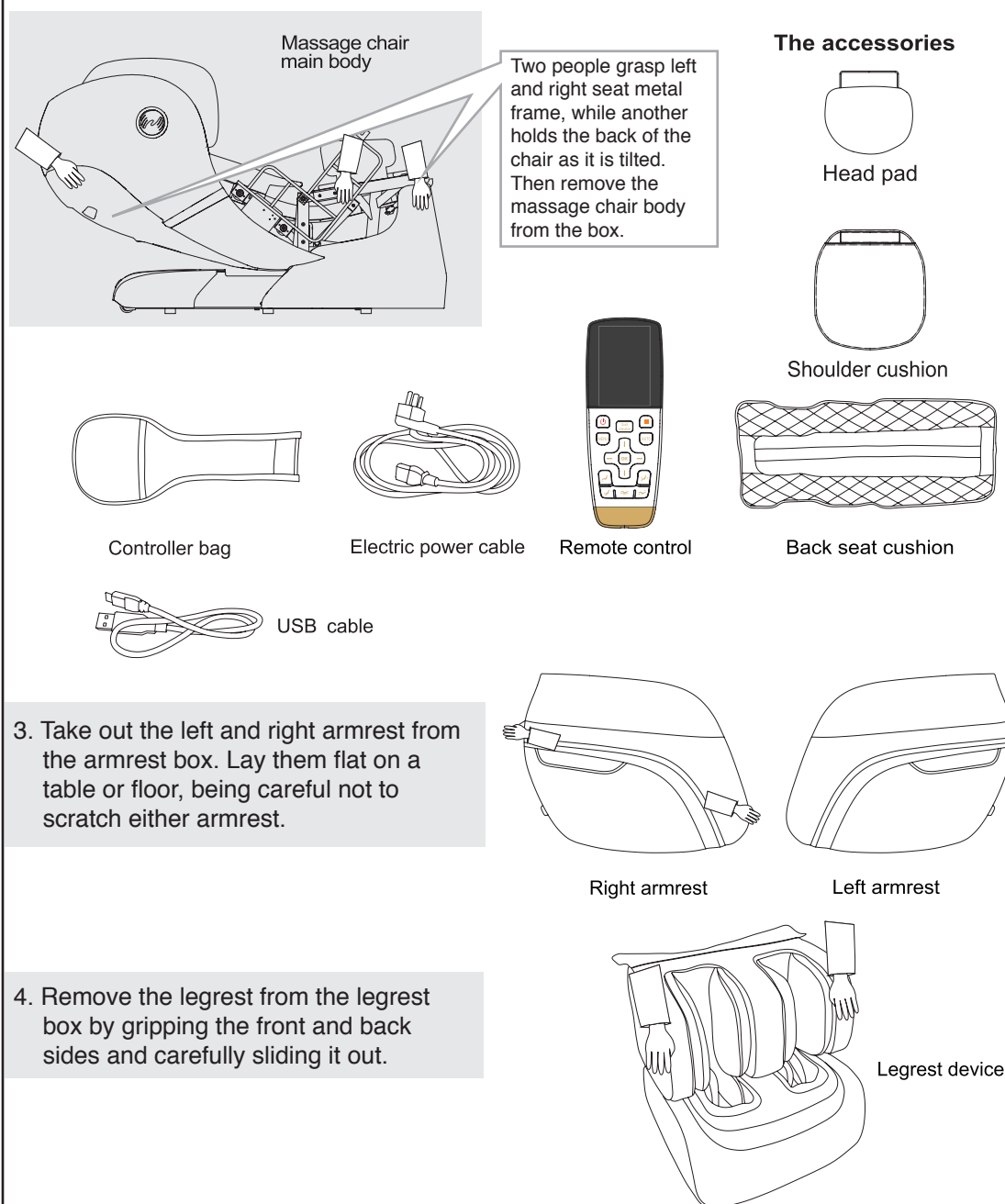


3. Take out the left and right armrest from the armrest box. Lay them flat on a table or floor, being careful not to scratch either armrest.

4. Remove the legrest from the legrest box by gripping the front and back sides and carefully sliding it out.

Installation (2)

When choosing a place for your new massage chair, check the surrounding environment and ensure there is enough space for the chair to recline fully.

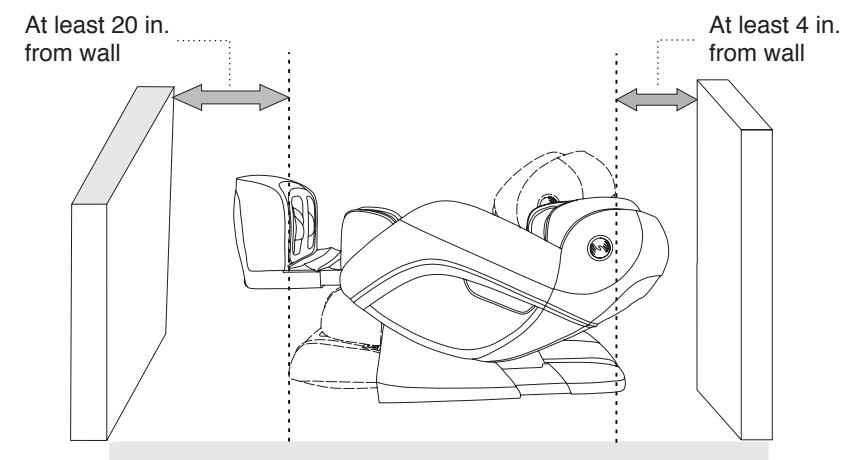

Please be careful when moving and using the chair, as it may scratch or damage flooring. We advise placing a small blanket, rug, or other protective covering to prevent damage to your floor.

Caution:


- Do not use the chair in wet places (e.g.: near a swimming pool, in a bathroom, etc.) to avoid risk of leakage or electric shock.
- Place the chair on a flat, level surface to avoid risk of falling or sliding and causing damage or injury.
- Do not place the chair under direct sunlight near heating appliances (e.g.: stove or oven) to avoid upholstery deterioration and damage.

Installation (3)

Step 2: How to assemble the legrest

1. Unscrew the bolt from legrest hanger.
2. Connect the legrest's airhose with the corresponding air nozzles of the main body tightly.
3. Connect power line to the power socket tightly.
4. Raise the rotation shaft board on hangers, lift the legrest, put rotation shaft into hangers on chair body, and lower rotation shaft board down.
5. Close the axis cover and screw on the bolt tightly.
6. Connect the leather of the legrest to the seat by securing the Velcro pieces.

CAUTION: Make sure the air tube and wires will not get tangled.

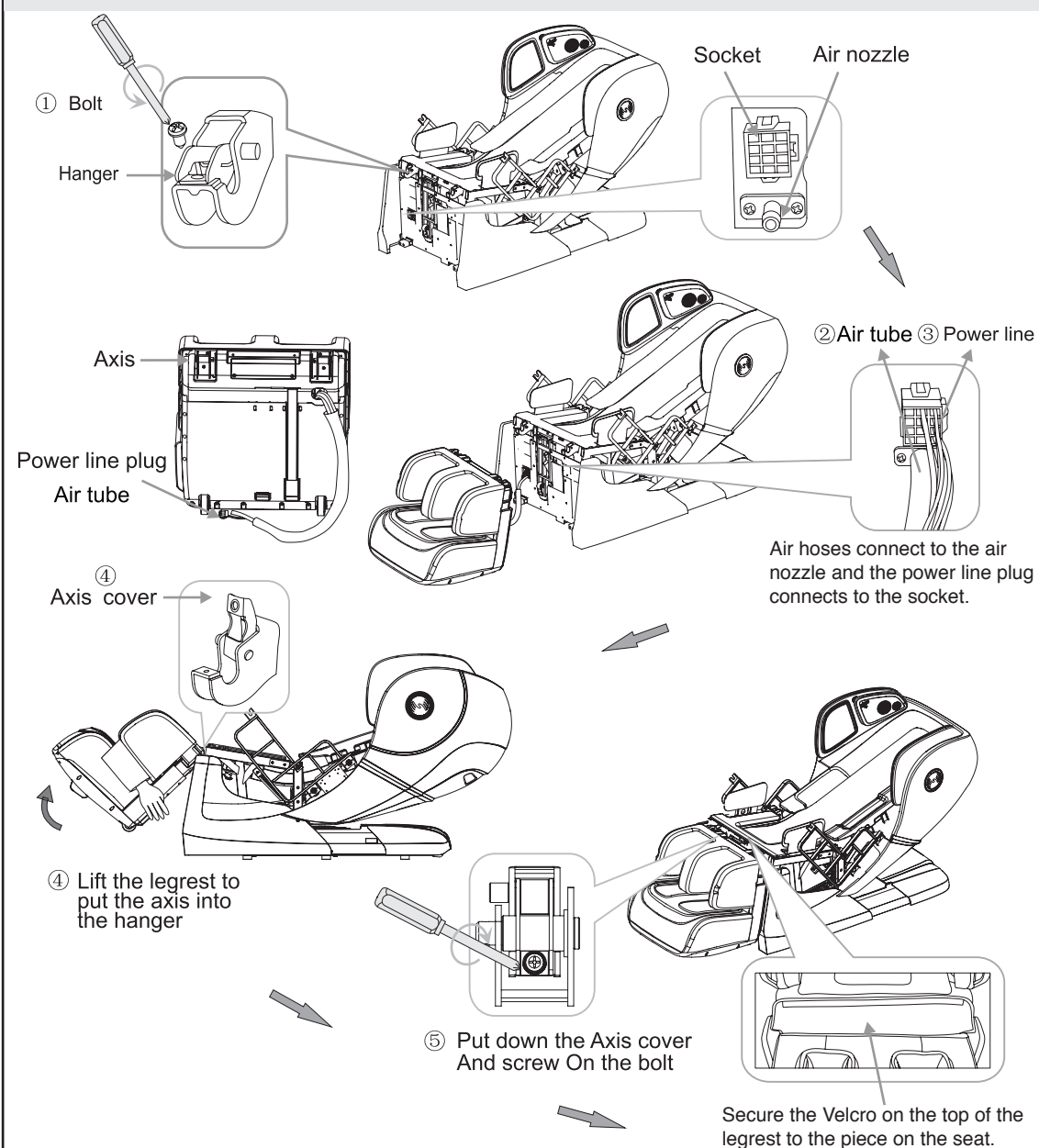


Installation (4)

Step 3: How to assemble the left and right armrests

1. Lift the armrest and plug it into the rotating stand on the side of the chair body.
2. Plug the power plug into the socket.
3. Connect the air hose of the armrest to the air tubes on the seat tightly. The black air tube should be connected to the air hose marked as "1" and the grey air tube should be connected to the air hose marked as "2".
4. Insert the rear armrest fixed plate into the fixing hole on the chair body.
5. Remove the rubber plug from the armrest, tighten the screws inside the hole, then put the rubber plug back on snugly.

CAUTION:




- Do not squeeze the air tubes during installation, as this may cause an air leakage or prevent the airbags from inflating.
- Ensure that the fixed plate and LED light wire are plugged in securely, otherwise the armrest may fall or the LED light will not work.

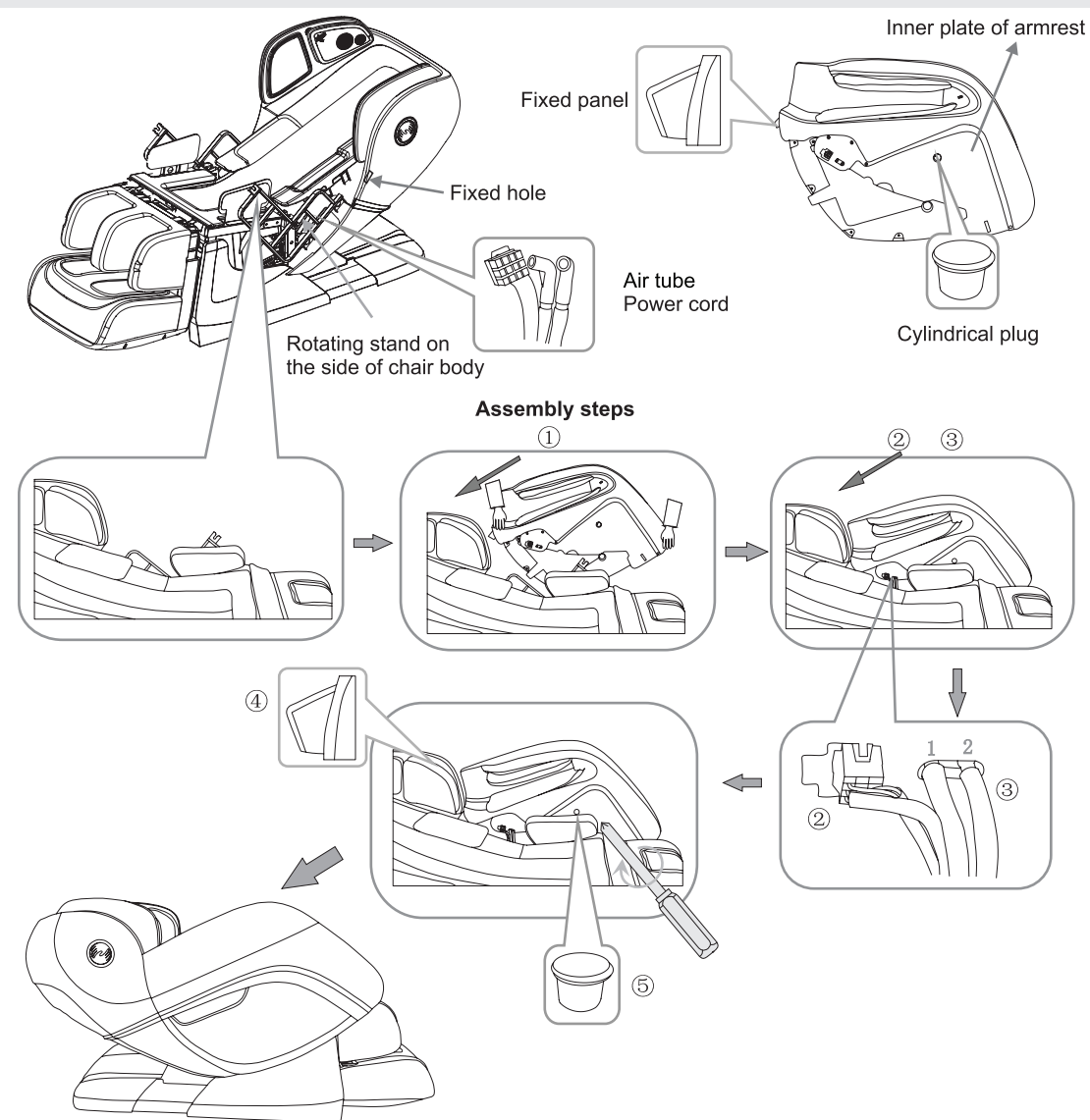

Installation (5)

Step 4: Install the back cushion, shoulder cushion, pillow, and remote control pocket

Put the back seat cushion, shoulder cushion, and head pad onto the headrest, then connect them with the zipper and Velcro.

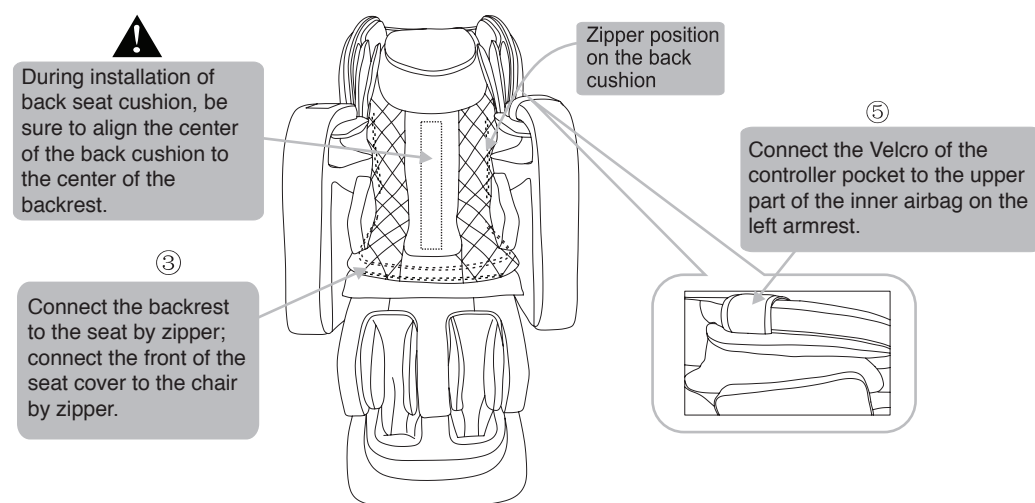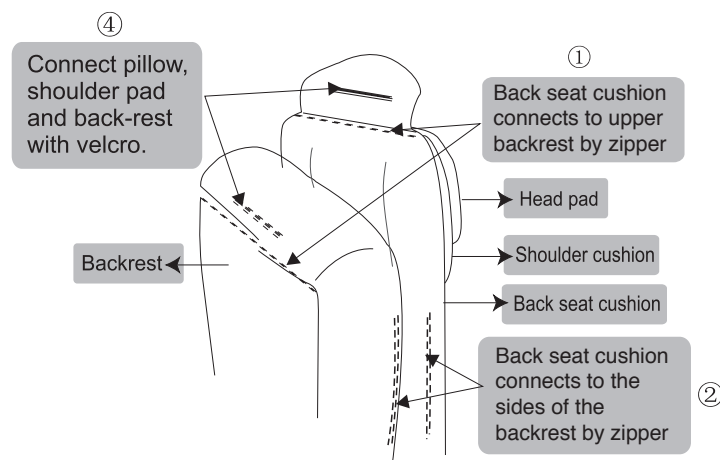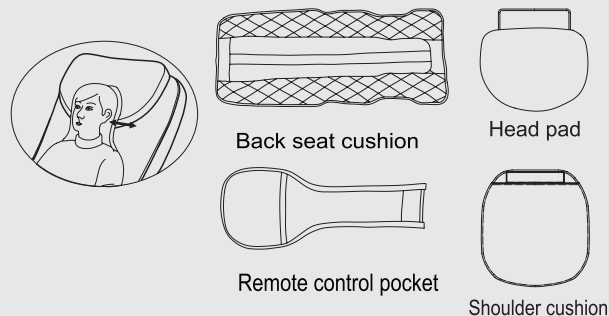


Installation (6)

Step 5: Restore the chair

1. Plug the power cord into the hole on the power box as shown in the right picture.
2. Insert the plug of the power cord into the outlet.
3. Turn on the switch behind the power box. Switch on/off "O/I". (Press "I" down to turn if on.)
4. Switch off, chair resets to upright state.
5. Use the convenient controller built into the armrest or the wireless controller to select your massage function.
6. When the chair is ON, the USB port can be used for charging your mobile device. This is also how you can charge the Overture's wireless remote, by using the USB cable that accompanies it.

⚠ Make sure the switch is OFF before plugging in the chair.

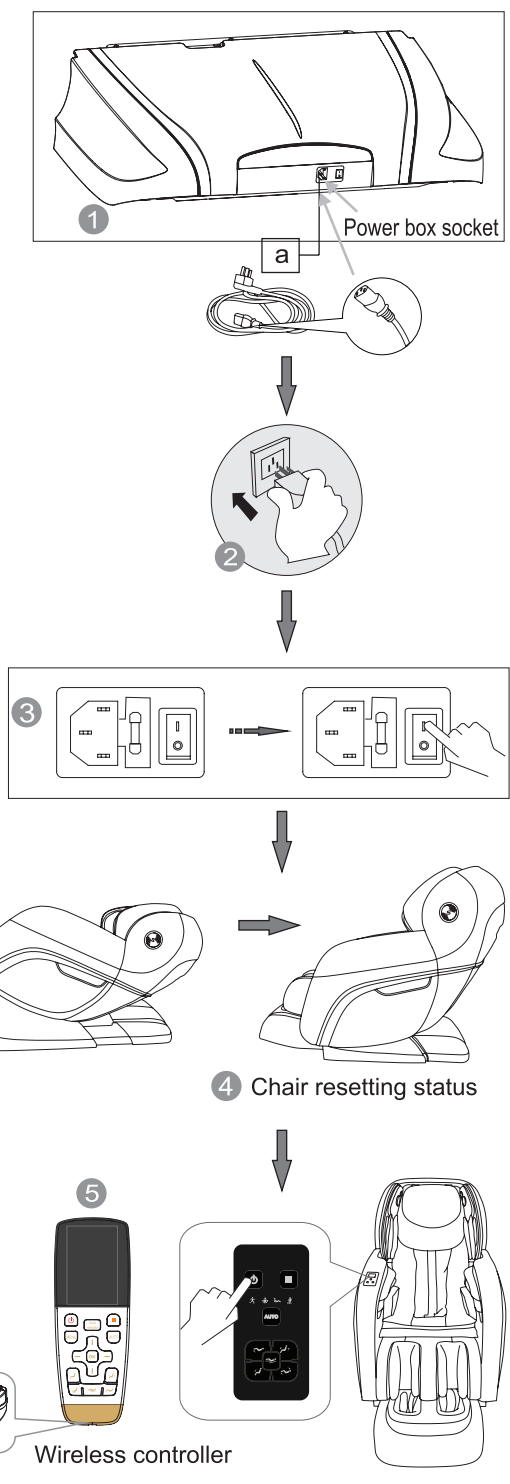

Environment



Caution



Don't place the chair in wet areas (e.g.: near a swimming pool, in a bathroom) to avoid leakage or electric shock.




Don't place the chair under direct sunlight or near a stove or other heating appliances to avoid upholstery deterioration or damage.



Place the chair on a flat, level surface to prevent the chair from falling or sliding and causing damage or injury.


● Check area around the chair

Ensure that there are no people, pets, or objects around the chair before reclining and extending the legrest.

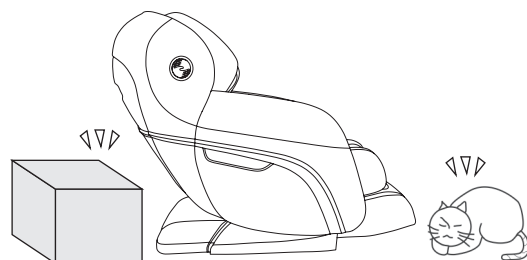

● Check power cable and plug


1. Clean the power cord and plug of any dust using a dry cloth.
2. Ensure the cord has not been damaged, bent, twisted, stretched, or tangled. Also make sure the chair is not on top of the power cord.

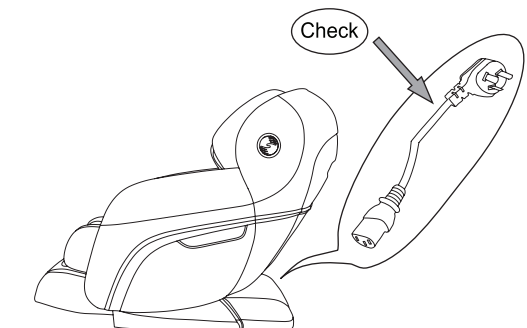



Check environment

● Where to place the chair

- At least 4 in. from the wall

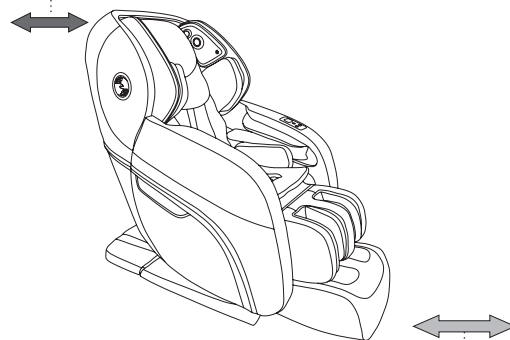

- At least 20 in. in front of the legrest

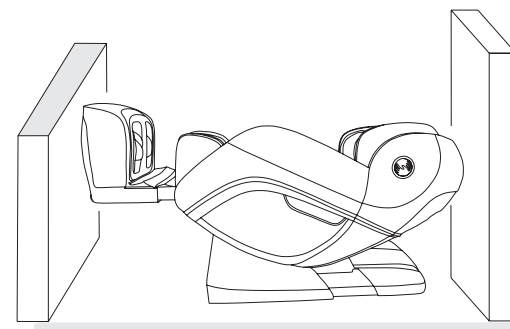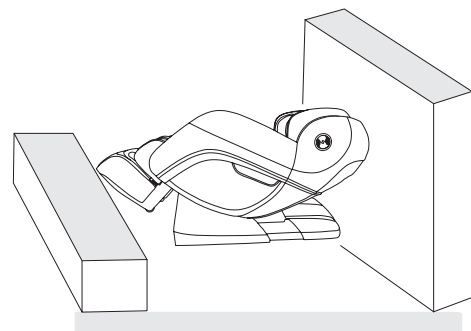

How to move massage chair



Caution




Don't move the chair during operation in order to prevent damage.

How to move the chair

- Unplug the cables and place them in the seat of the chair to avoid damage when moving.
- Do not slide on wooden floors or surfaces that are easily damaged. Use 2 people to lift and move the chair safely.
- Make sure the power is OFF and the plug is not dragging on the floor when moving.
- Keep the chair in the upright position when moving.



- 1 Press down on the head of the chair to tilt it back onto the castors. Then push the chair to the desired destination steadily and set it down slowly.

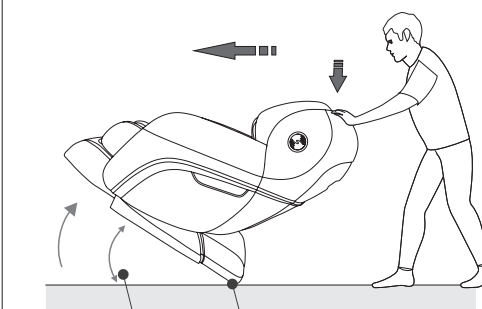

No more than 40 degrees — Rolling Wheel




Ensure the front of the chair is raised as shown in order to prevent any scratches on the footrest.

Please note: The armrest cannot be removed or adjusted when the chair is in the reclined position.

Adjustment before massage

Before using your massage chair:

- Make sure the massage rollers are in the starting position in the upper backrest, then sit in the chair slowly.
- During the body scanning function, if the massage roller mechanism cannot find the shoulder position, it will default to the middle of the backrest.

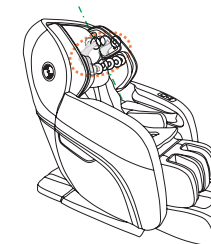


Correct position:

Sit up flush against the back of the chair, with your head on the pillow.




The remote control will beep once the body scan begins. When you hear the beep, you can press the up/down key on the control to manually adjust the shoulder position.

Try to remain still during the body scan in order to achieve an accurate and comfortable massage.

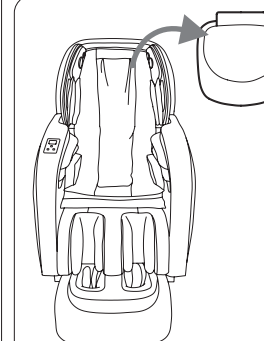


- If you feel the massage isn't strong enough, you can remove the pillow and shoulder cushion.

- We advise placing a small blanket, rug, or other protective covering to prevent damage to your floor.

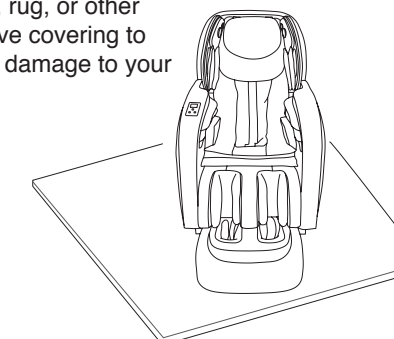

Built-in Remote

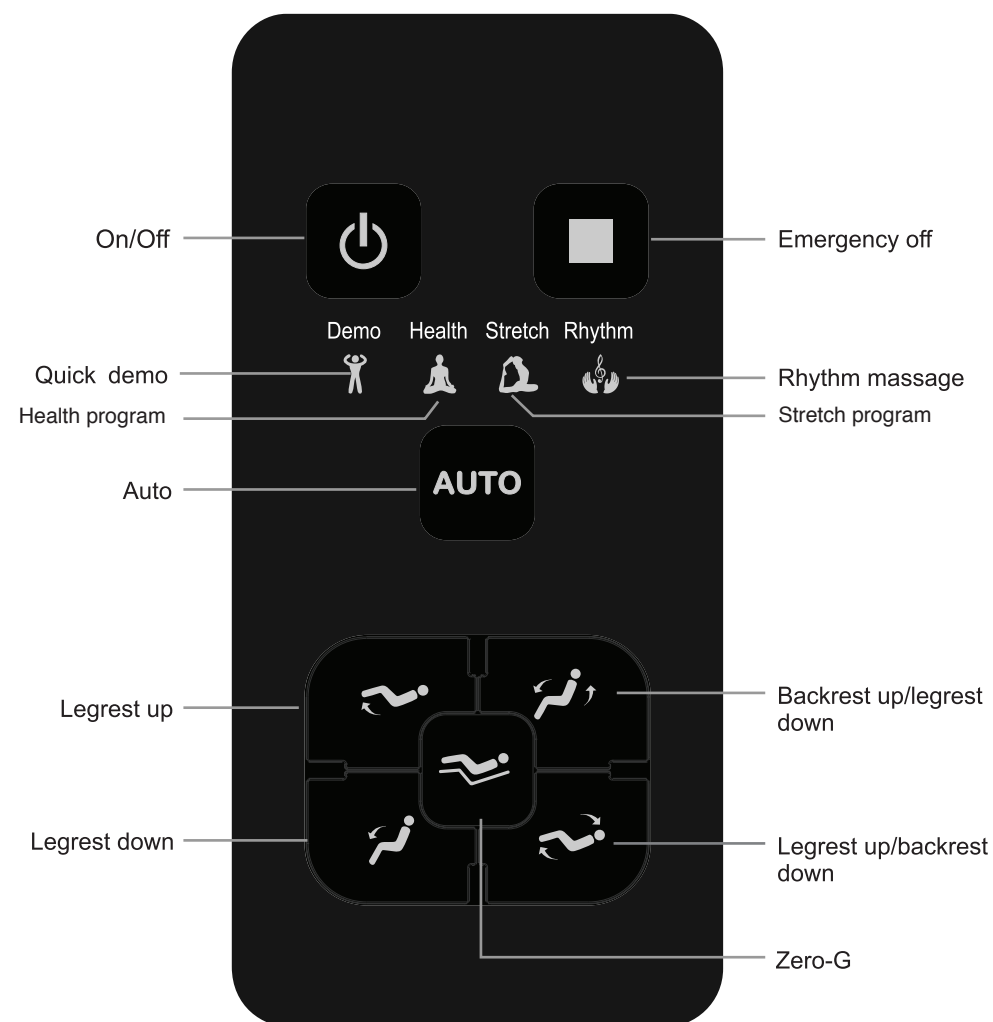

In standby mode, press the ON/OFF key to power on the chair. Choose a program through the menu on the wireless remote, or quickly select one by pressing the AUTO button on the built-in remote. Press the ON/OFF key again to turn off the massage chair, which will then reset to its upright position.



Press this key to stop all massage functions immediately. To start the massage again, you will need to press the ON/OFF key.



When the chair is powered on, use this key to toggle between four auto programs. The icon of each program will light up when selected.



Quick Demo: Demonstrates basic functions of the chair, default working time of 8 mins.



Health Program: Default working time is 20 mins.



Stretch Program: Default working time is 20 mins



Rhythm Massage: Default working time is 20 mins.



Hold down this key to raise the legrest; release the key and the legrest will stop moving. Please note: Legrest cannot be raised during "Stretch" program.



Hold down this key to lift the backrest and lower the legrest; release the key to stop. Please note: This function is not available during "Stretch" program.



Hold this key to lower the legrest; release the key to stop. Please note: Legrest cannot be lowered during "Stretch" program.




Hold this key to recline the backrest and raise the legrest; release to stop. Please note: This function is not available during "Stretch" program.



When the chair is ON, press this key once to set the chair at level 1 zero gravity; press again to reach level 2 zero gravity; and press a third time to reach sleeping mode. Press the key a fourth time to reset the chair to the upright position.

Wireless Remote

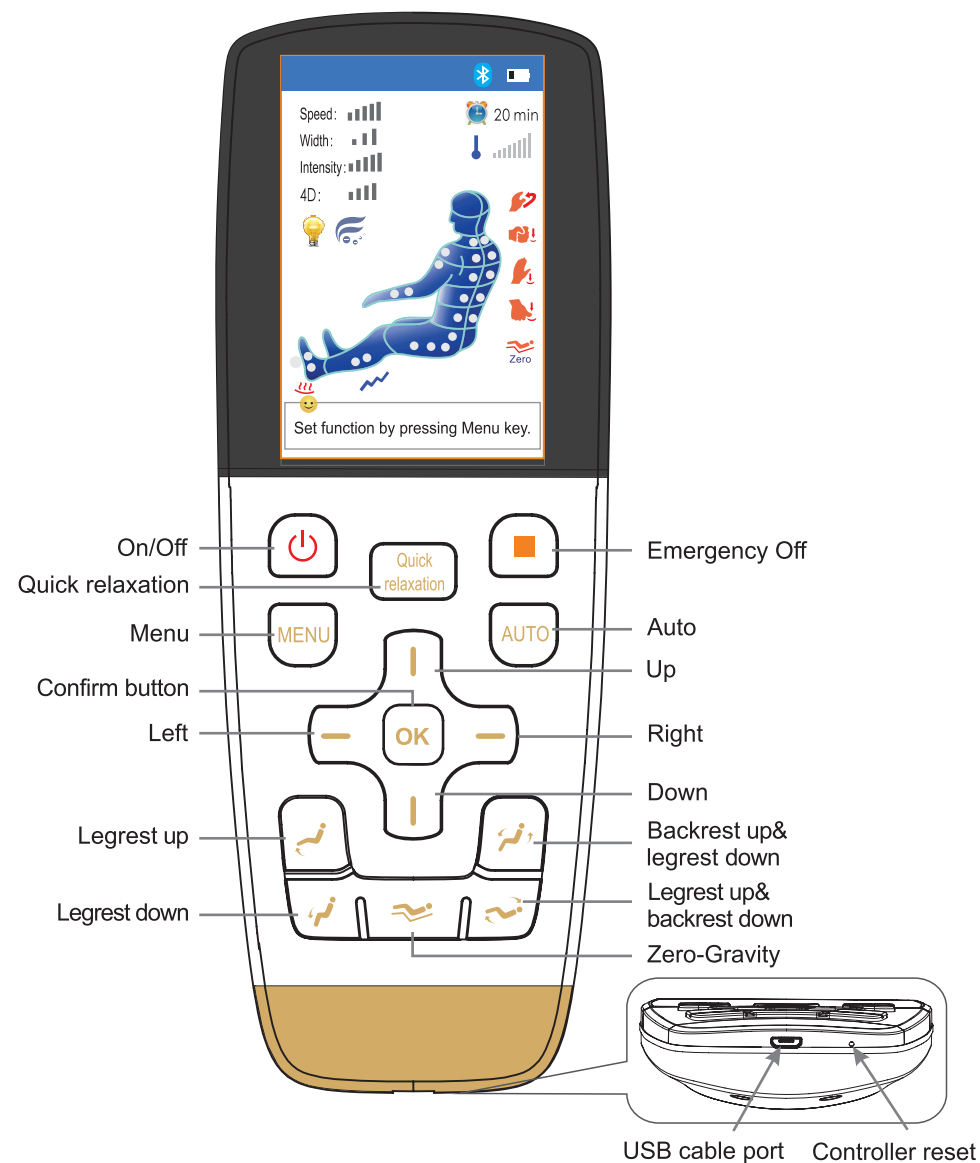**Connecting your wireless remote to Bluetooth for the first time**

1. Turn on the massage chair. Hold the MENU and OK buttons at the same time for 4 seconds, and the screen should indicate the Bluetooth pairing has begun.
2. The Bluetooth icon will light up if the connection is successful.
3. After pairing successfully, turn the controller off and on again. Your wireless remote is now ready to use to operate your massage chair.
4. The remote will remain paired and connect to your chair automatically every time you use it.

Please note: Your wireless remote should be charged upon receipt of your chair; however, if you find it won't turn on or has low battery, simply charge it with the provided USB cable.

Controller reset

If the remote disconnects or experiences pairing issues, use a pin or similar sharp object to press the reset button.



After switching the massage chair on, press the ON/OFF button for at least 4 seconds to fully turn on the chair. Then select from the AUTO programs or use the manual menu and keys to set a custom massage. Press the ON/OFF key again for 4 seconds to turn off the chair, which will then reset to its upright position. The remote will power off after about 3 to 5 seconds.



Press this key to enter Quick Relaxation mode. The massage status will display on the remote screen. Default working time for this program is 8 minutes.



Press this button during operation to stop all massage chair functions immediately. All the buttons will become invalid upon pressing the emergency stop button. Press the ON/OFF key to reactivate the remote control and re-start your massage.



Use this button when in the main interface to switch to the menu interface. Here, you can select different functions and settings for your massage. Press the Menu button again to return to the main interface.



Press this key to select an Auto program. There are six auto massage programs: Health, Relief, Stretch, Renew, Music Sync, and Calming. Your massage status will display at the top of the screen.



The direction key allows you to navigate up, down, left, and right. When in the Menu, use this key to select menu items. In Manual mode, use the up and down keys to adjust the position of the massage rollers during the shoulder height detection at the beginning of your massage and before a fixed point massage.



When in the Menu, use the OK button to confirm menu choices and enter the next menu. When on the main interface, press this button to turn on the air ionizer; press again to turn the ionizer off.



Hold down this key to raise the legrest; release the key and the legrest will stop moving. Controller screen will show "Position adjusting..." while adjusting the legrest.



Hold down this key to lift the backrest and lower the legrest; release the key to stop. Controller screen will show "Position adjusting..." while adjusting the massage chair.



Hold this key to lower the legrest; release the key to stop. Controller screen will show "Position adjusting..." while adjusting the legrest.




Hold this key to recline the backrest and raise the legrest; release to stop. Controller screen will show "Position adjusting..." while adjusting the massage chair.



When the chair is ON, press this key once to set the chair at level 1 zero gravity; press again to reach level 2 zero gravity; and press a third time to reach sleeping mode. Press the key a fourth time to reset the chair to the upright position.

LCD display graphic description

Main interface

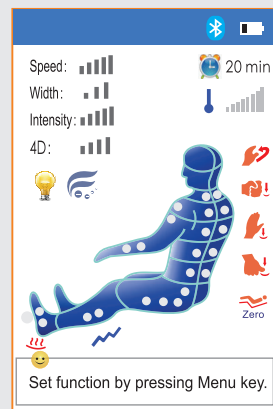

After switching on the chair, press the ON/OFF button and you will hear a short beep. The LCD screen will load and after about 2 to 3 seconds, it will display the main interface and enter standby mode. The main interface displays the massage position, method, time, current massage status, etc. All of this information will adjust once you start a massage.

Massage function

- Kneading
- Tapping
- Knocking
- Shiatsu
- Zero gravity
- Massage time left

Massage position

Massage Speed (Level 1, 2, 3, 4, 5)
 Massage Width (Narrow, Mid, Wide)
 Air Intensity (Level 1, 2, 3, 4, 5)
 4D Intensity (Level 1, 2, 3, 4)

Lumbar Heat Temperature Level


Heat On Indicator

Armrest LED light status

Air Ionizer Indicator


Leg rubbing function display

Menu interface

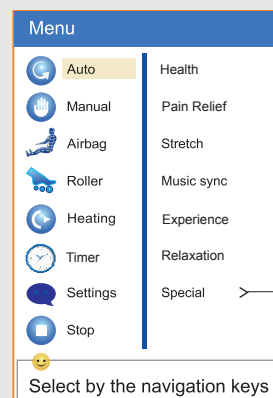


Press the Menu button when in the main interface to switch to the menu interface. Here, you can select all your massage functions and settings. Press the Menu button again to return to the main interface.

Auto Mode: Press Menu and select from Health, Relief, Stretch, Music Sync, Renew, and Calming. You can also select Special to view six additional programs: NeckShoulder, Waist & Spine, Back Rehab, Active, Relax, and Healing.

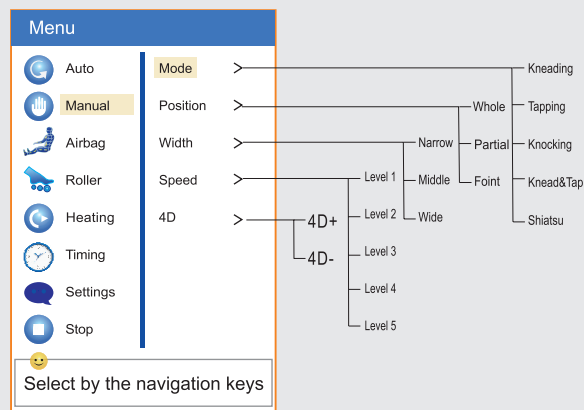

Manual Mode: Press Menu, select manual message using the navigation key. You have the option to adjust five settings: Mode, Position, Width, Speed, and 4D. Choose the corresponding function, then press OK to enter the next menu. Massage modes you can choose from include Kneading, Tapping, Knocking, Knead & Tap, and Shiatsu. Massage position lets you set the rollers at Whole, Partial (smaller area), or Point (precise spot).

LCD display graphic description

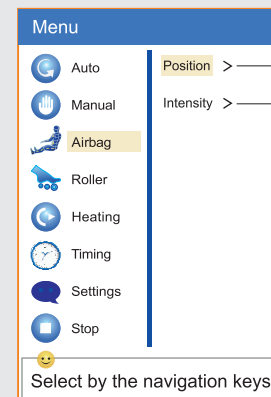

Airbag Massage: After pressing MENU key, navigate to Airbag, then press OK. You can adjust both the position and intensity. Simply choose the desired function, press OK, then choose your preference from the corresponding menu.

Foot Roller: Press MENU and navigate to Roller. Here you can choose from Quick, Soft, or Stop foot roller settings.




Heating: Press MENU, navigate to Heating, press OK, and select which heating functions you'd like to turn on or off. Options include: Back On, Foot On, Temp+, and Temp-.

Timer: Press MENU, navigate to Timing, and choose the length of your massage. Options include 10 minutes, 20 minutes, or 30 minutes. The default time is 20 minutes. The massage chair will turn off automatically and reset after the time is up.

Settings: Press MENU and navigate to Settings. Here, you can choose your language and turn the LED lights on or off.


Stop: Press MENU and select Stop to stop the current massage and return to the main interface.

When you begin a new auto program, the massage chair will begin the body scan and "Body scanning..." will display on the screen. After the body scan is finished, the chair will beep 5 times (once per second). During these 5 seconds, the massage rollers are in the shoulder position and you are able to press Up or Down to adjust. While adjusting, "Shoulder height adjusting..." will display on the screen.

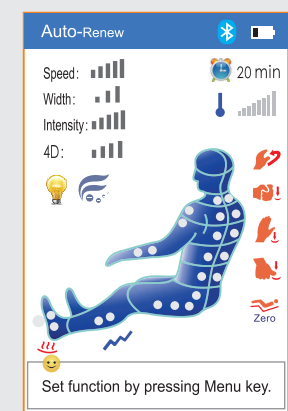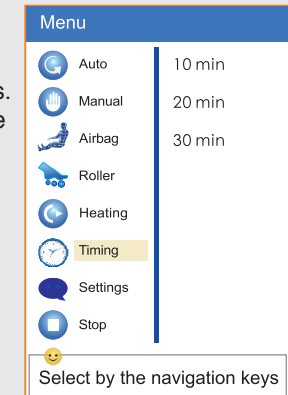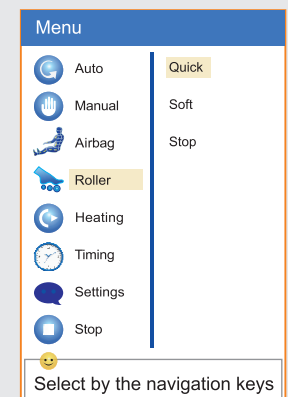

ON/OFF






Step 1: Power On

1. Plug the power cord into the port labeled "a".
2. Plug the cord into a power outlet.
3. Turn on the power switch on the power box.
4. Press ON/OFF button on the remote, which will beep once and then enter the main interface.

Step 2: Massage

1. Operate according to remote control operation instructions.
2. When the massage ends, the massage chair will return to the default position. During this time, "Shutting down" will display on the remote screen.

Convenient Built-in Remote Wireless Remote

Step 3: Turn Off the Power

1. Press the ON/OFF button to turn off the massage chair and it will reset automatically.
2. Turn off the power switch on the power box after the chair has fully reset and stopped moving.
3. Unplug the power cord from the outlet.

Function Illustration

● Legrest Up

Hold down this key to raise the legrest; release the key and the legrest will stop moving and stay in this position.




● Legrest Down

Hold this key to lower the legrest; release the key and the legrest will stop moving and stay in this position.



● Legrest up and backrest down

Hold this key to recline the backrest and raise the legrest; release to stop and keep the chair in this position.

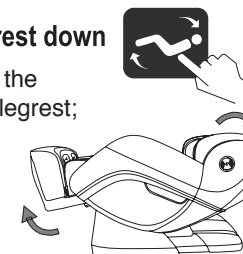


● Backrest up and Legrest down

Hold down this key to lift the backrest and lower the legrest; release the key to stop and keep the chair in this position.




● Legrest extend and retract

Legrest can extend up to 8.5 in. automatically.

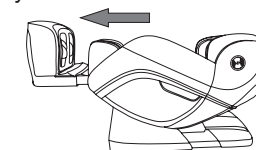


● Zero gravity

When the chair is ON, press this key once to set the chair at level 1 zero gravity; press again to reach level 2 zero gravity; and press a third time to reach sleeping mode. Press the key a fourth time to reset the chair to the upright position.

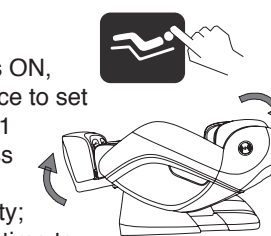


★ Zero gravity positioning was invented by NASA engineers. In this position, the chair reclines to elevate your knees above your heart, which allows the spine to decompress, provides a sense of weightlessness, and maximizes the effects of the massage rollers. Zero gravity position also helps increase blood circulation to and from the heart.

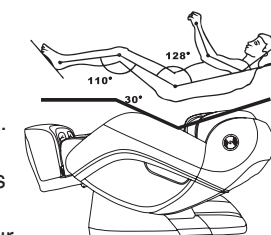

● New Curved Rail Provides Comprehensive Back & Buttocks Massage

The Overture features a 4D curved rail technology, which allows for both a highly effective decompression stretch and powerful glute massage. Unlike traditional L-track massage chairs, this roller track style allows the chair to recline to a flat position for a full stretch.

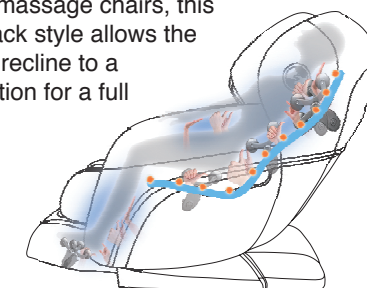


- When adjusting the backrest, beware of the gap between the backrest and armrest.
- Do not stand on or get up from the chair when the legrest is moving. Do not try to manually stop the legrest from moving, as this may cause damage to the chair.

Function Illustration


● Arm massage

The Overture contains capsules in the armrests with built-in airbags. Place your arms between the armrest's airbags to experience the kneading massage on your arms and hands.

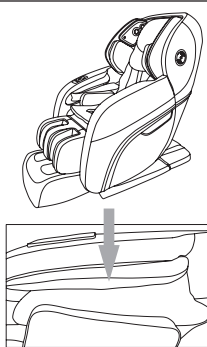


● Shoulder massage

Multi-layered airbags are located on both shoulders to provide kneading massage and to assist with the stretch function. The 4D massage rollers will push your body forward, while the airbags grip your shoulders to provide the ultimate full-body stretch.

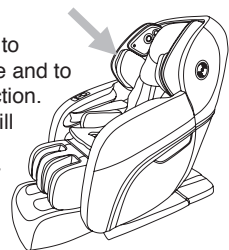

● Leg rubbing function

The Overture's innovative leg rubbing technology provides kneading and rubbing massage on the calves, which can help to relieve pain and rejuvenate the muscles.

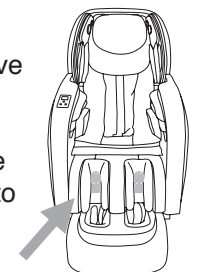


● Warm foot therapy

Powerful two-way foot rollers provide kneading massage to the soles of the feet. Combined with heat therapy, these foot units provide a truly soothing massage to tired, hard-worked feet.





● Ionizer Function

The chair's air ionizer uses air purification technology, which captures the harmful substances in the air and purifies the air quality. This feature can help improve sleep quality, boost metabolism, and strengthen the immune system.

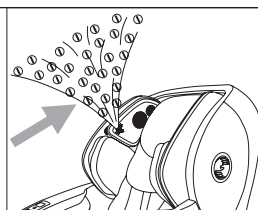

● Music therapy

Relax to your favorite sounds with the Bluetooth-enabled speakers located in the Overture's headrest. You can even sync your massage to the beat of your music with the Music sync program.

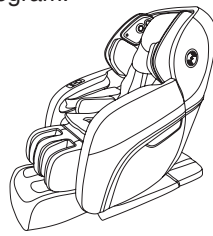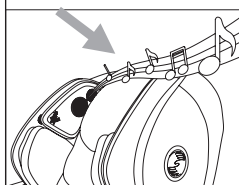

● Voice Operation

To use the Overture's speech recognition feature, first use the command, "Turn on speech recognition mode." It will then be able to identify the rest of your commands. You will be able to hear the corresponding speaker voice feedback after the chair successfully recognizes your command. The corresponding responses are as follows:


COMMANDS	RESPONSES
Turn on speech recognition mode	Speech recognition has been activated
Voice off	Voice recognition mode is off
Massage chair shut down	Massage chair is powering off
Health massage	Health massage has been activated
Pain relief	Pain relief massage has been activated
Body stretch	Body stretch massage has been activated
Zero gravity start	Zero gravity position
Experience	Experience massage has been activated
Relaxation	Relaxation massage has been activated

How to Use Voice Operation:


1. Say, "Turn on speech recognition mode." If you say "Voice Off" and "Massage chair shutdown," the massage chair will shut down by default. To restart, say "Turn on speech recognition mode" again.
2. When your voice command is recognized, the speaker will respond after 1 second. You can then say your next command entry.
3. After turning on speech recognition mode, the chair will remain in this mode for 3 minutes. If you say nothing after 3 minutes, it will shut off speech recognition mode and you will need to say "Turn on speech recognition mode" again.
4. Use standard language and speak slowly and clearly.

Restoration of massage chair


- 1 Make sure there are no people, pets, or objects within the range of the backrest and legrest.

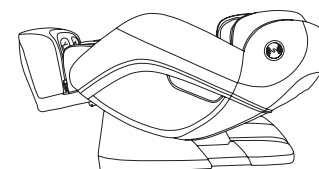


- 2 Restore the backrest and legrest to the normal position using the Backrest Up & Legrest Down" button, or press the On/OFF key to rest the chair to its default position.

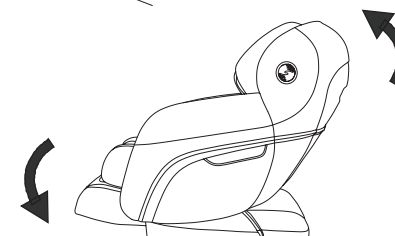


- 3 After the chair returns to its upright position, store the wireless remote in the pocket on the left armrest.

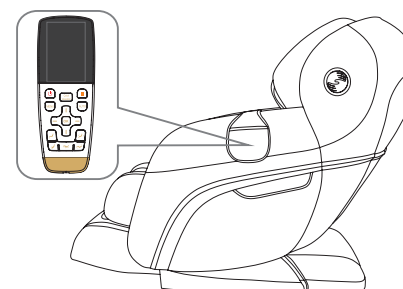

- 4 When getting out of the chair, hold the armrests with two hands and stand up slowly, moving your feet out of the legrest. Then step onto the ground. Restore the backpad and pillow to their original position.

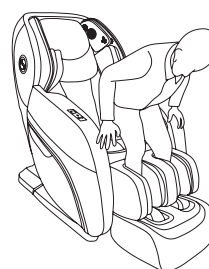


✓ Right

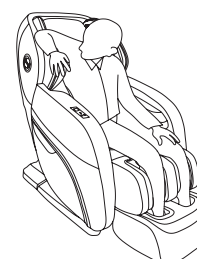


✗ Wrong

- 5 Switch off the power.

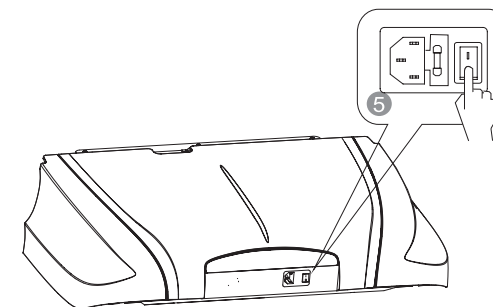

- 6 Unplug the power cord from the outlet. Make sure to grip the plug properly, as shown below.

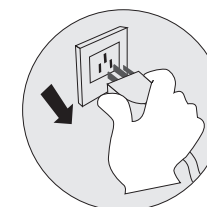

- 7 Clean the chair if necessary.



Attention!

- After using the chair, press ON/OFF to end all massage functions and wait for restoration.
- Make sure to cut the electricity to avoid any unexpected injuries.

Synthetic leather

● General cleaning

Clean with a soft, dry cloth.



Do not use chemical cleaning products, as they could damage the chair upholstery.



Forbidden



● If leather is dirty

① Dampen a soft cloth into 3-5% neutral detergent solution of gentle soap and water. Wring it out, then clean the stain on the surface gently.

② Wipe with a water-dampened cloth.



③ Wipe the area with a dry cloth.

④ Allow to air dry. Do not use a blow dryer.

Plastic parts

① Dampen a soft cloth into weak neutral detergent solution of gentle soap and water. Wring it out, then clean the stain on the surface gently.



② Wipe with a water-dampened cloth.

③ Allow to air dry.


Clean the controller

① Wipe with a water-dampened cloth.

② Allow to air dry. Do not use chemical cleaning products..



Forbidden

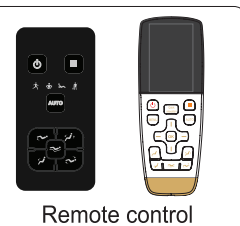

Remote control

Cloth material

① Dampen a soft cloth into weak neutral detergent solution of gentle soap and water. Wring it out, then clean the stain on the surface gently. Do not use chemical cleaning products.

② Gently brush the dirty spot(s) with the weak neutral detergent. Avoid scrubbing excessively.

③ Wipe with a water-dampened cloth.

④ Allow to air dry.



Cleaning Precaution

Before cleaning, unplug the electricity. Don't touch the plug with wet hands in order to avoid electrical shock or fire risk.



FAQ and Troubleshooting

Trouble	Possible reason	Troubleshooting
Massage chair won't work even when powered on.	Check if the control is on or not.	Turn on control switch
	No massage function has been chosen	Choose massage function
	Plug is not securely in the outlet.	Make sure plug is securely connected to the outlet.
	Wire or plug is damaged.	Contact Infinity for repair and replacement by an authorized service provider.
	The fuse burned out.	Replace fuse of same specification.
	Something is wrong with the inner circuit.	Contact Infinity for repair and replacement by an authorized service provider.
Massage chair is making a beeping sound while operating.	It is the air pump, motor, and mechanical parts.	This is a normal sound.
Left and right massage rollers are at different heights.	Massage rollers alternate during operation.	This is a normal function.
Massage chair suddenly makes loud noises during operation.	This could be a result of using the chair for a long period of time.	Turn off the power and let the massage chair rest for at least a half hour before using again.
	Inner components have become worn down due to long-time, heavy use.	Contact Infinity for repair by an authorized service provider.
Massage chair suddenly stops during operation.	The chair has been powered down suddenly.	Reconnect to power and turn it back on.
	Setting time is up.	Turn off the power and let the massage chair rest for at least a half hour before using again.
Neither the backrest nor legrest will recline or lift/extend.	Check for obstacles around the chair. If the chair overloads, it will shut down automatically.	Remove the obstacles; turn off the power and allow the chair to rest for a half hour before using again.
The chair won't return to its original position.	Excessive operation; setting time is up.	Contact Infinity for repair by an authorized service provider.
Electrical wire or plug is emitting heat.	It may be overloaded.	Turn off the power and let the massage chair rest for at least a half hour before using again.
	Other abnormal case	Stop use and contact Infinity for repair by an authorized service provider.



Please Note: If the issues you're experiencing are not solved by any of the above troubleshooting options, please turn off the massage chair completely and unplug it. Then contact Infinity Massage Chairs directly for assistance and repair by an authorized service provider.

Technical data

Name	INFINITY	Power	280W
Model	Overture	NW/GW	Main body: 253/302 lbs Armrest: 32/43 lbs Legrest & Footrest: 66/76 lbs
Input voltage	<div><div></div>AC120V, 60Hz<div></div>AC220V, 50/60Hz<div></div>AC240V, 50Hz</div>		
Safety design	Class I	Packing size	Main body: 59.8 x 31 x 38" Armrest: 38.5 x 15.5 x2 8" Legrest & Footrest: 23 x 20 x 24" (L x W x H)