




*The World's Best Massage Chairs*



---

68 Route 125  
Kingston, New Hampshire 03848  
U.S.A.  
T. (603) 347-6006 F.(603) 642-9291  
[www.infinitetherapeutics.com](http://www.infinitetherapeutics.com)

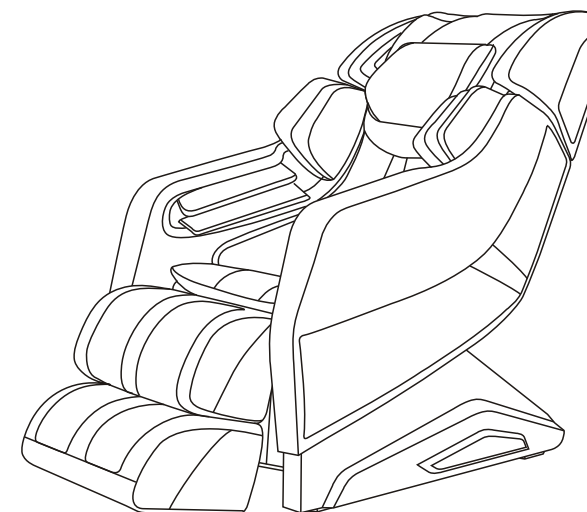

**Deluxe Multi-functional Massage Chair**

**IT-Riage-X3  
USER MANUAL**

## ► TECHNICAL SPECIFICATIONS

---

Name: Deluxe Multi-functional Massage Chair

Model: IT-Riage-X3

Rated voltage: 100-120V~

Rated frequency: 50/60 Hz

Rated power: 240W

Safety design: Class I

Thank you for your purchase of this massage chair. Before use, please read this USER MANUAL carefully, and pay special attention to the safety precaution to ensure proper operation and use.

Please keep this USER MANUAL properly for further reference.

## ► CONTENTS





---

SAFETY PRECAUTIONS .....	1
NAME AND FUNCTION OF EACH COMPONENT .....	3
External structure .....	3
Internal structure .....	3
Power box .....	4
Controller .....	5
Controller buttons .....	6
LCD screen icon decomposition .....	7
FUNCTIONS AND FRATURES .....	12
PRECAUTION BEFORE USING .....	13
Installation Method .....	13
Move method .....	17
Turn on power supply .....	17
Installation position .....	18
Floor protection .....	18
Massage chair commissioning .....	19
Before sitting in the chair .....	19
Blue tooth .....	20
After using .....	20
Surrounding environment .....	20
Grounding .....	20
PRODUCT STORAGE AND MAINTENANCE .....	21
TROUBLE SHOOTING .....	22
TECHNICAL SPECIFICATION .....	23




## SAFETY PRECAUTIONS

### Safety Precautions

- ▶ Before using, please read this “SAFETY PRECAUTIONS” carefully to ensure proper operation.
- ▶ The following “SAFETY PRECAUTIONS” is of great help to your correctly use of this massage chair.
- ▶ To clearly indicate the degrees of danger and damage, the classes of accidents caused by incorrect use are classified into “WARNING” and “CAUTION”, instructions marked as followings should be strictly followed:

	<b>WARNING</b>	Improper operation may cause injury to users.
	<b>CAUTION</b>	Improper operation may cause users injury or product damage.
	Forbidden. Specific instruction will be indicated in the diagram drawings.(E.g, the left drawings means disassembly is not allowed.)	
	Compulsory Specific instruction will be indicated in the diagram drawings.(e.g, the left drawings means to disconnect the power plug)	

### WARNING

 COMPULSARY	We recommend that those with any of the following conditions should consult a doctor before using.
	1. Those who are using electronic medical instrument imbedded in the body, such as cardiac pacemaker etc.
	2. Those that are being treated by doctors, especially those who feel uncomfortable.
	3. Those who are malignant tumor patients, heart disease patients, acute patients.
	4. Women who are pregnant or menstruating.
	5. Those with osteoporosis or fracture spine.
	6. Those with skin disease or those whose skins are injured.
 PROHIBITED	7. Those whose body temperature is over 38 degrees (Febrile Phase).
	Do not allow the mentally disabled to use the chair alone. (It is OK under supervision)
	Do not allow children play on massage chair.
	Do not allow more than one person to use this massage chair at a time.
	If the power cord or power plug becomes damaged in any way you must disconnect from the mains supply immediately and consult an authorized service engineer for repair or replacement.
 NO DISASSEMBLY	Do not damage, process, excessively bend, pull hard, twist or knot power cord.
	Juveniles are not allowed to use this product.
	Anyone else, except maintenance personnel specified by our company, should not disassemble, maintain, or remodel this massage chair. If any problem, please contact local dealer immediately.

## TROUBLE SHOOTING

No.	Problem	Possible Cause	Solution
1	The sound “Xu Xu”, “Ta Ta” is heard.	It is a normal sound made by the working motor or the sound is produced when the massage ball is rubbed against the cloth cover.	It is normal and no need for treatment.
2	The machine can not be started up.	The power supply plug is detached.	Insert the power supply plug into the socket again.
		The power supply switch is not turned on.	Turn on the switch
		The fuse of massage chair blew out.	Replace a new one with the same mode
3	The airbag installed at seat cushion or legrest can not work.	To check whether there is air tube blockage.	To clear the air tube

If abnormal phenomenon still remains after above measures have been taken to the massage chair, please contact local distributor for repair and maintenance.

**Note:** Do not dismantle the product by yourself, and our company shall shoulder no responsibility to the damage of the product or the injury of people caused by improper use! Thanks your cooperation.

## ▶ PRODUCT STORAGE AND MAINTENANCE

### Maintenance

- ▶ Please twist the cotton cloth which has been dipped in neutral detergent before wiping the back cover, hand rest, legrest, air tube and some other plastic components, then use a dry cotton cloth to wipe them again.
- ▶ Use dry cloth to wipe the controller, DVD, and power supply box.
- ▶ Use a piece of cloth that has been dipped with neutral detergent to wipe the back cushion and pillow, then use a piece of wet cloth to wipe them, and finally put them to a well-ventilated place for natural drying.

#### ⚠ Caution

- ▶ Make sure to pull out the power supply plug firstly before maintenance. Do not use wet hands to pull out the plug to avoid electric shock. When clean the back cover, please be cautious of detergent entering into the power supply box!
- ▶ Do not iron the upholstery.

### Storage method

- ▶ Wipe off dust and dirt from the product after use, if the product is not to be used for a long period of time, please cover the product with a dust cover (self provided).

#### ⚠ Caution

- ▶ Do not store the product under direct sunlight or in places of high temperature. This may cause the color of the upholstery to fade.

### Kind reminding

During daily use, please check whether the following symptoms exist:




- ▶ Scorched smell exists.
- ▶ When touch the wire, the power is on, but sometimes off.
- ▶ The power supply cord is hot.
- ▶ Other abnormal symptom

#### ⚠ Caution

- ▶ Please stop using the massage chair to avoid breakdowns and accidents occurring if any above symptom exists.
- ▶ Turn off the power supply switch and pull out the power supply plug, and then contact the local distributors or headquarters for repair and maintenance.
- ▶ Except all the normal operation and maintenance mentioned above, any other repair and maintenance can only be operated by our authorized dealer or service agent.

## ▶ SAFETY PRECAUTIONS


### ⚠ CAUTION

 COMPULSARY	<ul style="list-style-type: none"> <li>◆ Please do not continuously use the massage chair over 20 minutes, and it is limited to massage a particular part of body in 5 minutes.</li> <li>◆ Please put the massage chair on flat ground for use.</li> <li>◆ When adjusting the legrest or backrest, please do not exert excessive force on it; in addition, keep children and pets well clear.</li> <li>◆ When removing power plug, please hold the insulation plug, do not pull power cord.</li> <li>◆ Please do not use the power supply that is beyond the specified electric voltage, as this may cause fire disaster or damage massage chair.</li> <li>◆ Before using, please uncover the back cushion to check whether backrest or fabric is damaged; if damaged, please stop using this massage chair immediately, and contact franchiser for repair immediately.</li> <li>◆ Please do not operate the controller with damp hands, and do not put heavy objects on controller.</li> <li>◆ Only people less than 120kgs are acceptable.</li> <li>◆ Please do not use other therapeutic equipments together with this massage chair.</li> <li>◆ During use, if you feel uncomfortable, please stop use immediately and consult your doctor.</li> <li>◆ This massage chair is designed and intended for non commercial household use only, do not use the massage chair for medical treatment.</li> <li>◆ The massage chair is equipped with heating function, so people who are insensitive to heat must use it carefully.</li> <li>◆ Always switch off power and remove the plug from the mains socket immediately after use and before cleaning.</li> </ul>
 PROHIBITED	<ul style="list-style-type: none"> <li>◆ To avoid causing discomfort, do not use the massage chair immediately after meals.</li> <li>◆ Do not fall into sleep on massage chair, when the massage chair is working.</li> <li>◆ Do not use the massage chair when you are drunk or feel bad.</li> <li>◆ If the upholstery becomes dirty, please use only known good quality leather cleaner, do not use chemical cleaners, petrol or detergents etc, and do not spray pesticide on it.</li> <li>◆ Do not sit on backrest, legrest, or armrest to avoid accidental injury or machine fault.</li> <li>◆ Please do not smoke while using the massage chair.</li> </ul>
 REMOVE PLUG	<ul style="list-style-type: none"> <li>◆ Before maintenance, please pull out power plug. In addition, please do not insert or unplug the power plug with damp hands, to avoid electric shock.</li> <li>◆ When the massage chair is not in use, please switch off power and remove the plug from socket. When power is off, please cut off power supply immediately and pull out the power plug from socket.</li> </ul>




## ► NAME AND FUNCTION OF EACH COMPONENT

### External structure

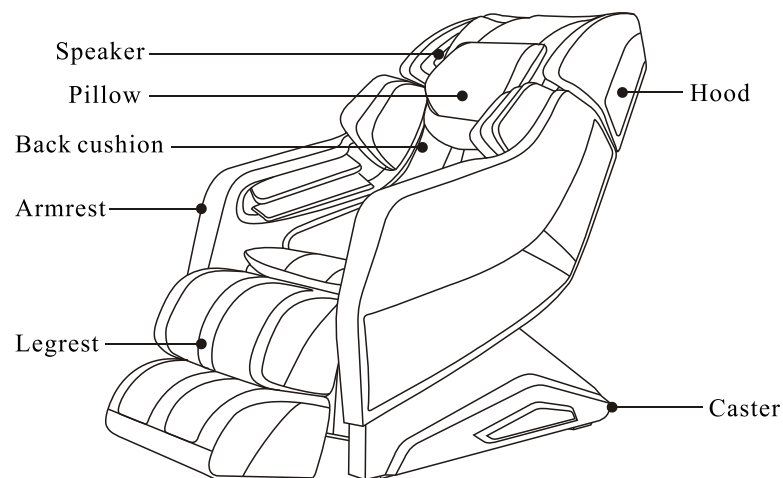

### Internal structure

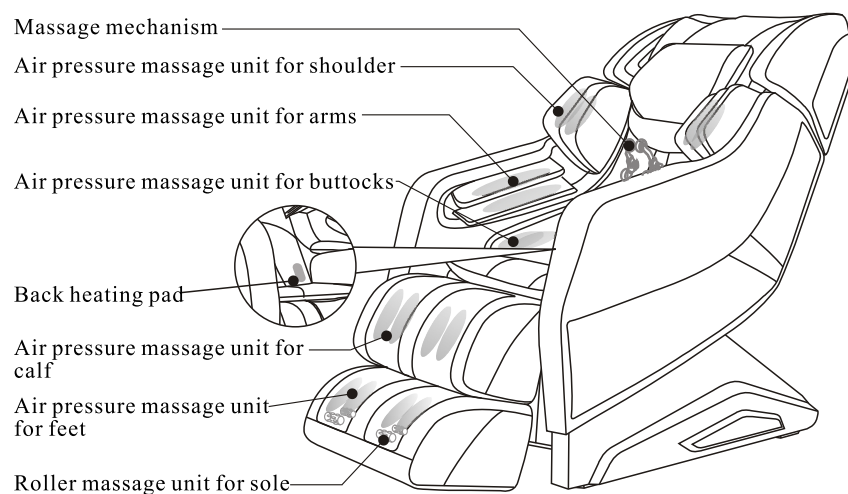

## ► PREPARATION BEFORE USING

### Blue tooth

- This product has blue tooth connecting function via audio devices. Open blue tooth from audio devices, match well with IT-Riage-X3 for music play.

### After using

- Please power off the chair and unplug the chair power cable from the power supply socket after using well the massage chair.
- Please store the chair if it is not in use for a long time.

### Surrounding environment



Keep the massage chair away from damp places (for example: bathroom) to avoid electric leakage.



Keep the massage chair away from high temperature sources (for example: stove) to avoid fire or damage to leather material.


### Grounding

This massage chair belongs to Level One electrical appliance, make sure to use three-core power supply socket and to guarantee the ground wire connected to the power supply socket has obtained good ground connection to avoid electric leakage, electric shock and some other negative effect during use.

## ► PREPARATION BEFORE USING

### Massage chair commissioning

- Before commissioning, please make sure there is enough space around the chair. Then use the controller to recline down the chair backrest to the limit position, rise up the legrest to the highest position and stretch out to the limit length.
- Testing all functions work normally, then power off the chair, return back to original position.  
(How to select the program, please refer to “Controller instruction” and “display screen menu instruction”)

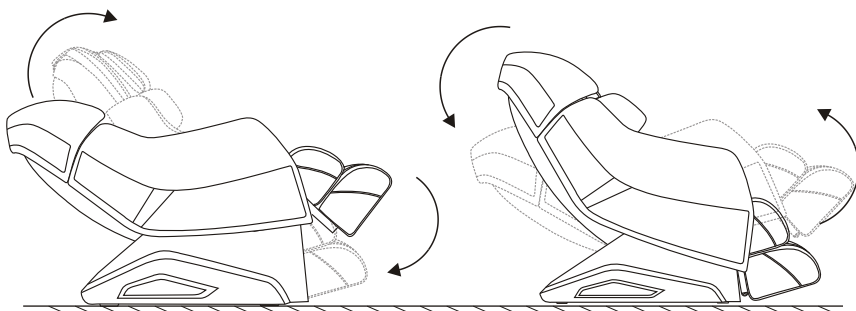

#### ⚠ Caution

- Do not put fingers or any other objects into the gap between the seat and legrest.
- When adjusting up and down angles of the backrest and legrest, please make sure there is no obstacle around the massage chair.
- During using, please make sure the power cord is not entangled or clamped but laid orderly.

### Before sitting in the chair


- Make sure the legrest is well installed.
- Make sure the chair is in original position, if any part is not in initial position, please power on/off again to make the whole chair go back to original position.

#### ⚠ Caution

- If the legrest rises up, please do not seat on the chair.
- Do not make the legrest to suffer from heavy loading suddenly, or it may damage the legrest mechanism.

## ► NAME AND FUNCTION OF EACH COMPONENT

### Power box

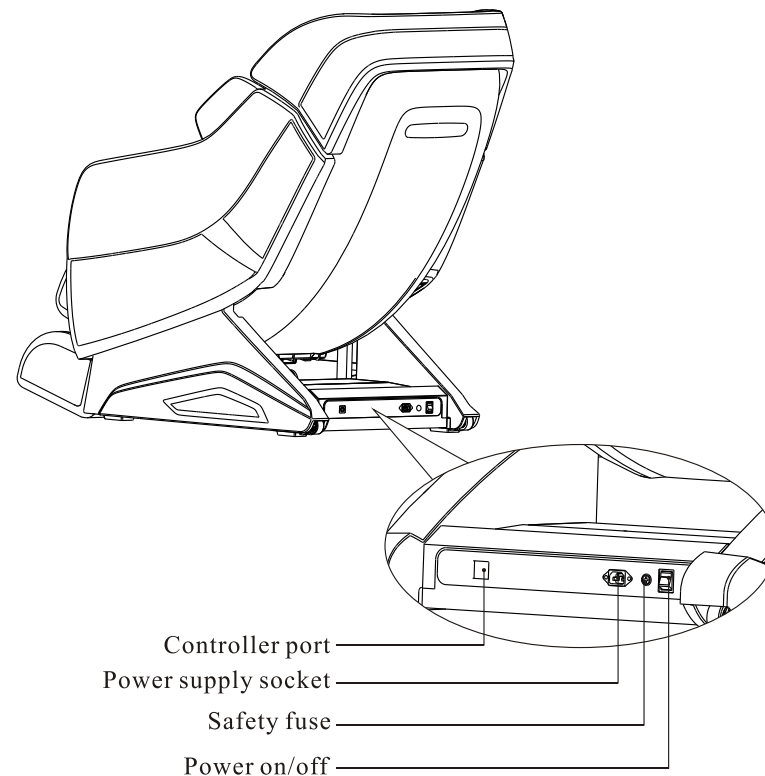


#### ⚠ Caution

- Do not let children play on the chair in case that accident happens due to the power supply box.
- Please power off the chair and unplug the chair from the power supply socket after the massage finishes
- Do not use the chair under thunders in rainy days.
- Please power off the chair and unplug the chair from the power supply socket when the chair is not in use, so that children won't accidentally get the chair started.
- Please store the chair if it is not in use for a long time.

## NAME AND FUNCTION OF EACH COMPONENT

### Controller

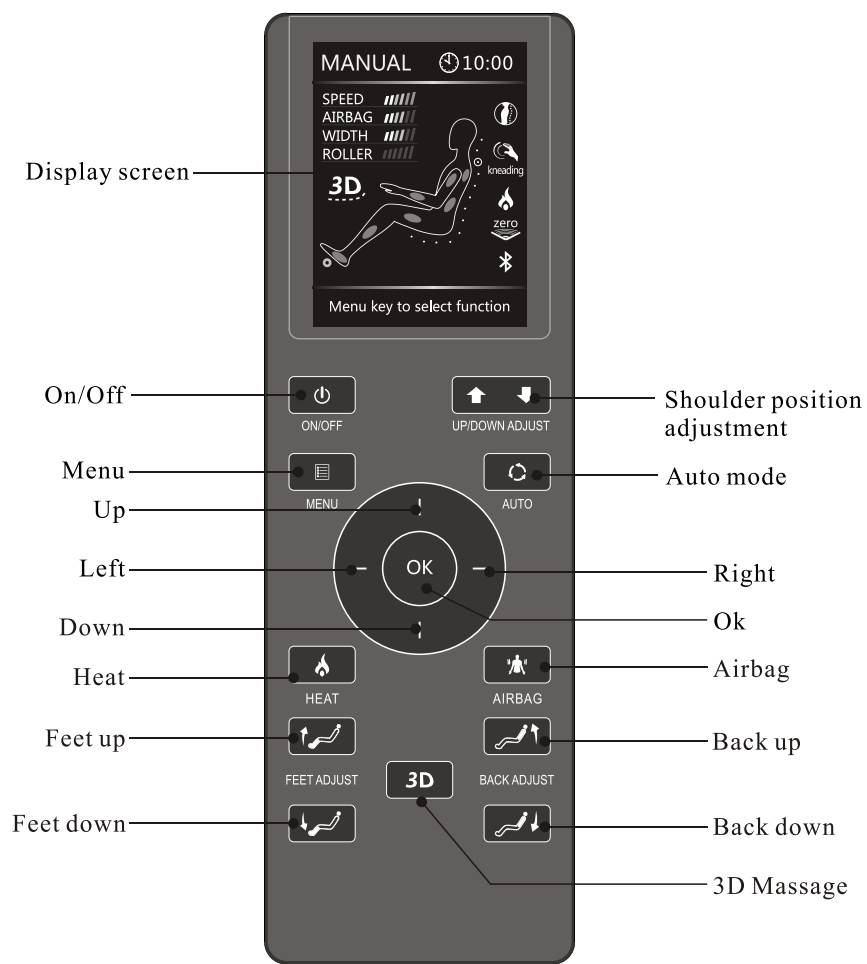


## PREPARATION BEFORE USING

### Installation position


Safety distance for the periphery of massage chair.

- Please make sure there is adequate space for installation this product:  
Back: minimum 60cm from the wall or any other blocks;  
Front: minimum 30cm from any blocks;
- In order to prevent from signal interference, please keep at least 1m away from the TV, radio and other video and audio devices.

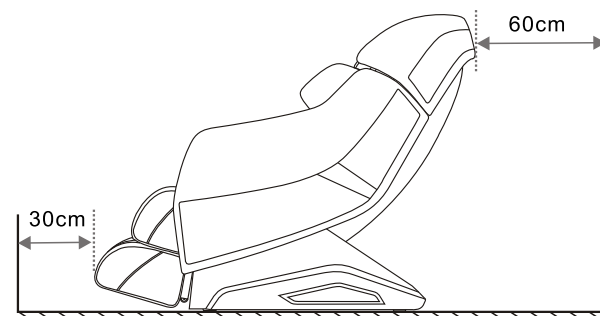

### Floor protection


- Please put carpet or other pad under the massage chair, in order to avoid any unpredictable impression on the floor.

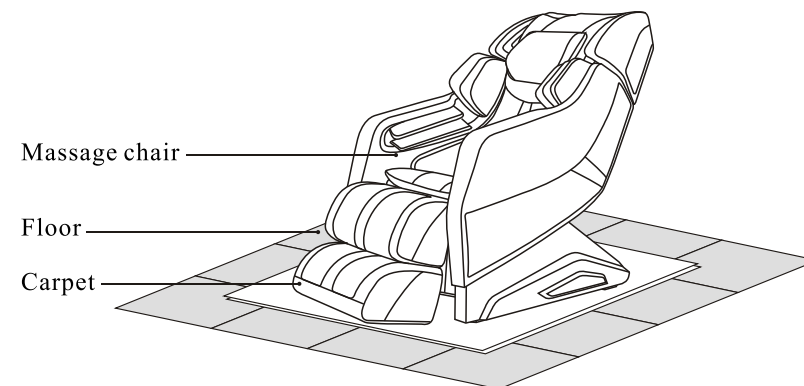

## ► PREPARATION BEFORE USING

### Move method

- There are moving rollers at the bottom of the chair.
- Raise the back to the highest point.
- Switch off and plug out the chair.
- Lift the front part of the chair as showed in the figure(one person raise the front part via holding moving holder, another person push the chair when holding the top of the chair) , move the chair by the rollers.

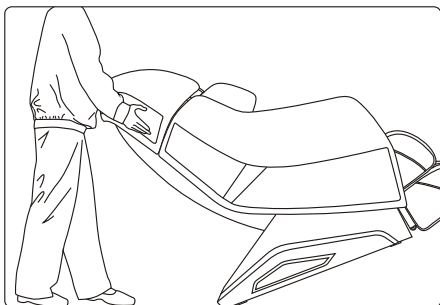


### ⚠ Caution

- Don't move it when there is people in it.
- Don't use the armrest while moving the chair, only hold the moving holder in the back and front

### Turn on power supply

- Plug the power cord into three core socket with well grounding.
- Tun on the power switch on the bottom of massage chair.

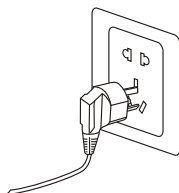

### ⚠ Caution

- Before switching on the chair, please make sure there is no damage on the power wire and other exposed wires.
- Make sure the power switch is off before the chair is connected to electricity.

## ► NAME AND FUNCTION OF EACH COMPONENT

### Instruction of Controller

#### ⏻ On/Off

- When power is on, click this button, display screen is lightening, then click up/down button to choose the auto mode, manual mode, air pressure massage, heat, 3D massage, and set menu.
- Click this button again to turn off all functions. The default working time is 20 minutes

#### ☰ Menu

- Click this button to choose auto mode, manual mode, air pressure massage, heat, 3D massage, and set menu.



#### Shoulder position adjustment

- Under manual and point mode, press UP/DOWN to adjust the position of point massage.



#### Auto programs

- When massage chair is working, press auto button to enter auto program menu, you can choose sports refresh, extention, rest& sleep, working relief, neck&shoulder or waist&spine according to your own body state and massage time.(3D massage excluded). (See the LCD menu for detail)



#### Direction adjustment

- Under menu and auto mode, press the button to adjust direction. When program is determined, press the middle OK button to confirm it and enjoy the chosen massage.
- During the massage, press OK button to switch over zero gravity, Press up or down button to adjust mechanism speed. Press left or right button to adjust airbag force.

#### 3D 3D massage

- Under manual mode, press 3D button to adjust the 3D force.



#### Heating





- Press this button to start the heating function,the indicating in the screen will be lighten;Press this button again,heating function will be shut off,the indicator in the screen will be off.



#### Airbag

- Press this button to enter the adjusting program for air pressure massage,press up and down button to adjust massage part and air pressure strength.


## NAME AND FUNCTION OF EACH COMPONENT

- 
**Feet up**  
 ▶ Press this button, the legrest will be moved up to the top position, release this button, it will stop.
- 
**Feet down**  
 ▶ Press this button, the legrest will move down to the lowest position, release this button, it will stop.
- 
**Back up**  
 ▶ Press this button, the backrest will be moved up to the original position, release this button, it will stop.
- 
**Back down**  
 ▶ Press this button, the Backrest will be moved down to the lowest position, release this button, it will stop.

## Instruction of Display screen menu

### ▶ Auto mode

When electricity is connected, press power button, LCD is on. Press auto to enter auto program menu, you can choose sports refresh, extension, rest& sleep, working relief, neck&shoulder or waist&spine according to your own body state and massage time. (See the menu picture below )

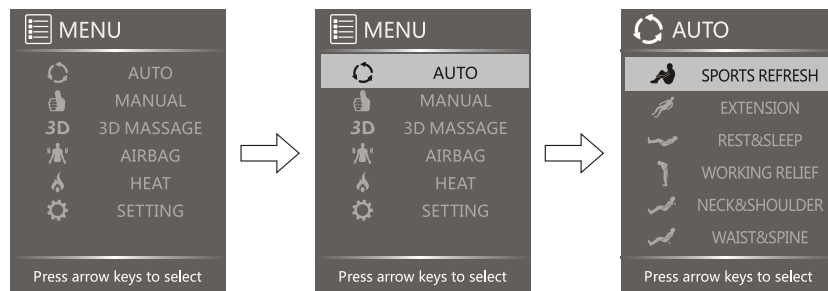

**Sports refresh:** This is designed for the muscles after sports activities and exercise. It can ease and relax muscles after intense exercise, especially promote the body circulation system to improve the muscle recovery after exercise.

**Extention:** This massage comes from Thai massage. the strong kneading and stretch of the body can effectively stretch legs, waist and other parts, and it is the best choice to relax tired muscles and increase physical vigor.


**Rest and sleep:** To improve the quality of lunch break and sleep; The mild massage kneads body parts from heavy strength to light one, from light to no strength to relieve tired body gradually.

**Working relief:** This is an exclusive program for white collars, who work long time at desk and often make business trip, it can ease muscles and recuperate body.

**Neck and Shoulder:** To relieve neck and shoulder pain symptom. 80% of the massage time will focus on the neck and shoulder. The curve design of the chair backrest makes the massage rollers completely touch the shoulder and neck. So it can effectively relieve your fatigue of these parts.


## PREPARATION BEFORE USING

- ▶ Lift the legrest unit, and connect the connectors of legrest to the relevant ports on the front frame of the seat frame well. (see diagram below)
- ▶ Connect the Pu covers of the legrest unit to the seat with zipper.

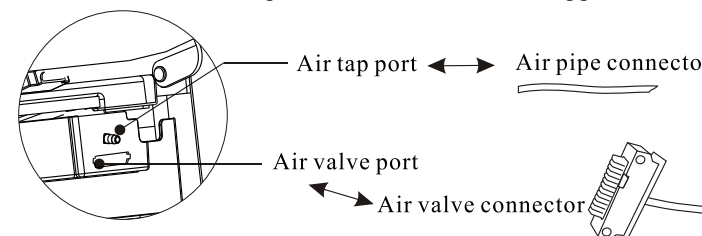

### Step 4 Installation of other units

- ▶ Insert the plug of remote control into the socket on the power box at the bottom.
- ▶ Connect the power cord connector to the power cord port on the power box. (see diagram below)

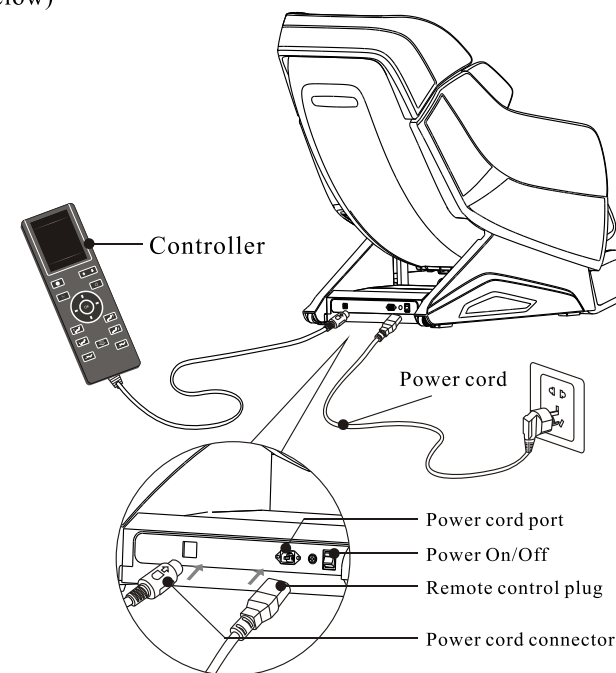


- ▶ After all parts assembled, please insert the power cord into the power socket, and then turn on the switch, press the power on the remote control to start massage. (see diagram above)


## ► PREPARATION BEFORE USING

### Step 3 Installation of legrest unit

- Lift the calf part and connect the calf with the seat frame by hanging (pay attention to the balance), move the hanger to the axis on one side, then do the same for the other side.

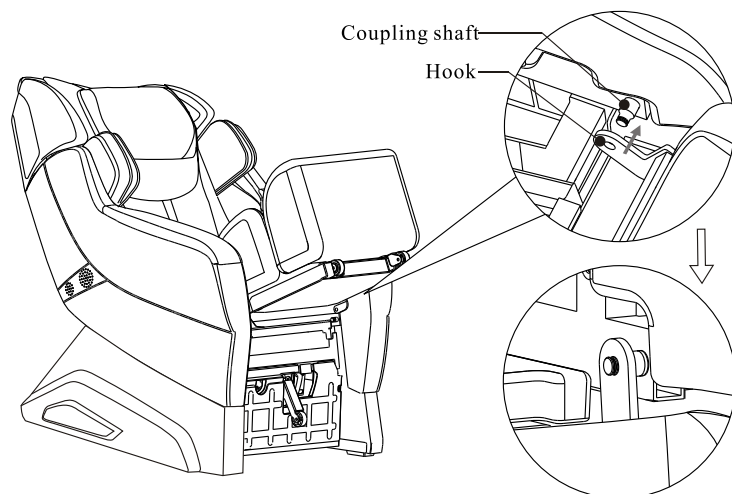

- After the calf is well connected with the axis, fasten it by the lock on the left and right side. ( see the figure below)

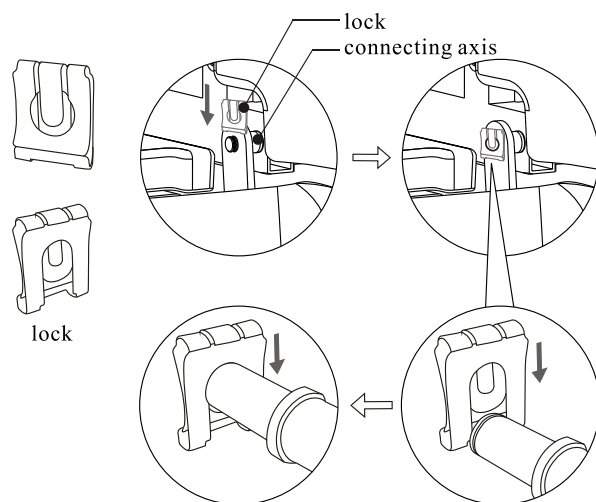


## ► NAME AND FUNCTION OF EACH COMPONENT

**Waist and Spine:** To relieve the pain of lumbar vertebra. 80% of the massage time will focus on the lumbar vertebra and the rollers will conduct specific massage to relieve the sense of tension and keep the health of this part.


### ► Manual Programs:

Press Menu button, choose Manual and press it to enter manual programs. Press up or down button to choose mode, part, width, speed, roller or zero.

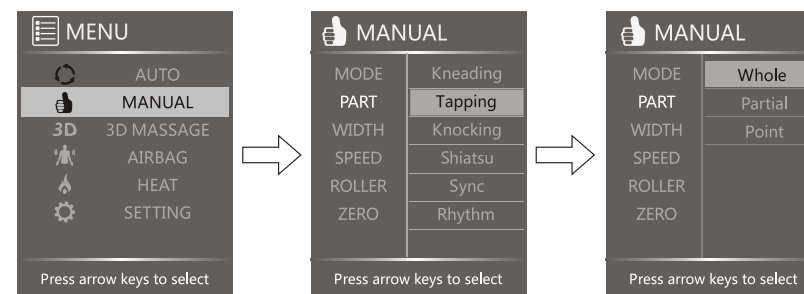

**Mode:** Under mode function, press right button first, then press up or down button to circulation choose massage manipulation of kneading, tapping, knocking, shiatsu, sync, or rhythm.

**Part:** Under Part function, press right button first, then press up or down button to circulation choose whole, partial or point. Back massage mechanism is by whole, partial and point.

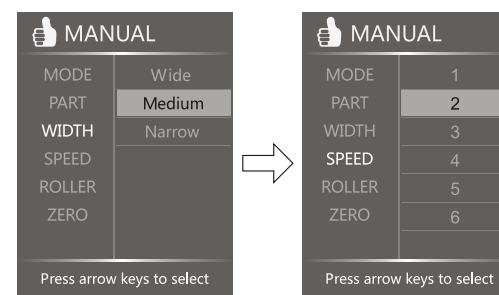

**WIDTH:** Under Width function, press right button first, then press up or down button to circulation choose Wide, Medium or Narrow.


**SPEED:** Under SPEED function, press right button first, then press up or down button to circularly choose speed from one to six.

**ROLLERS:** Press MENU button first, then choose ROLLER, adjust the roller speed from 1 to 3 or stop the rollers by pressing up or down button.


**ZERO GRAVITY:** Zero gravity creates a completely relaxing, naturally comfortable massage position: Somatological design makes back, buttocks and leg stay in a certain angle like lying in space capsule, under which the body pressure is thoroughly released and the muscle come into completely relaxing status.



## NAME AND FUNCTION OF EACH COMPONENT

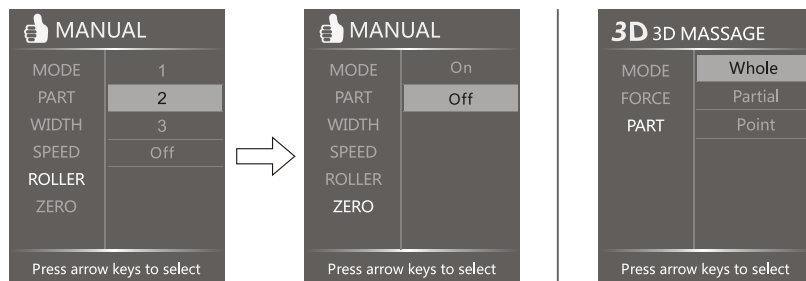

- **3D MESSAGE:** Press MENU button first, then choose 3D, adjust the 3D MODE (Deep shiatsu, healthy breath, massage extend), massage force and part.

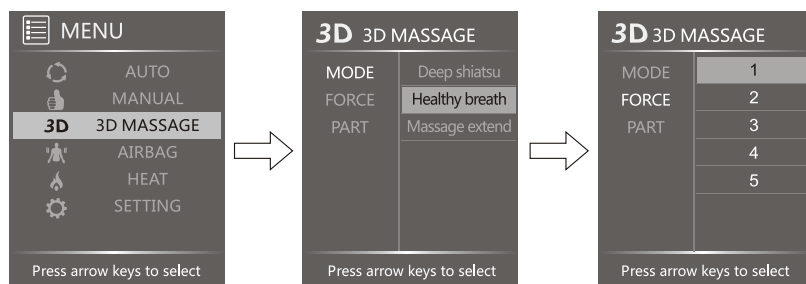


**Deep shiatsu:** deep shiatsu on the back for 10 minutes;  
**Healthy breath:** 5 minutes breath health works along with chair mechanism's rhythm for breathing in and out, which does good to the user's health.  
**Massage extend:** 5 minutes massage extend.

**FORCE:** Under manual, 3D massage force is adjustable from 1 to 5 grades.


**Part:** Under Part function, press right button first, then press up or down button to circulation choose whole, partial or point massage.

- **AIRBAG:** Press MENU button first, then choose AIRBAG, press up or down button to choose airbag place and force.

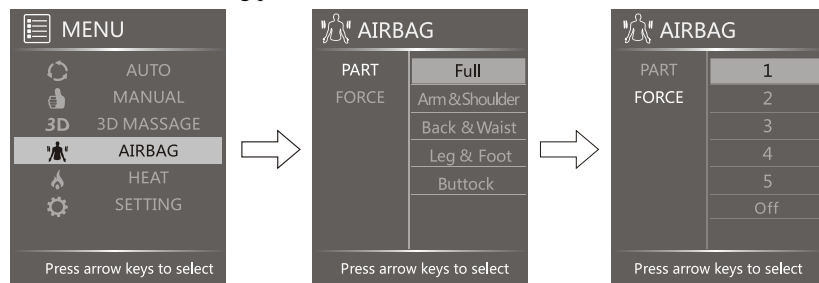

## PREPARATION BEFORE USING

- Hang the hook of the armrest well to the connecting rod of the backrest and then connect the armrest to the side of the seat frame, push it down to make it connected well. (Diagram 2)
- Fasten the armrest to the seat frame with one screw on each side. (Same assembly method for each side of armrest)

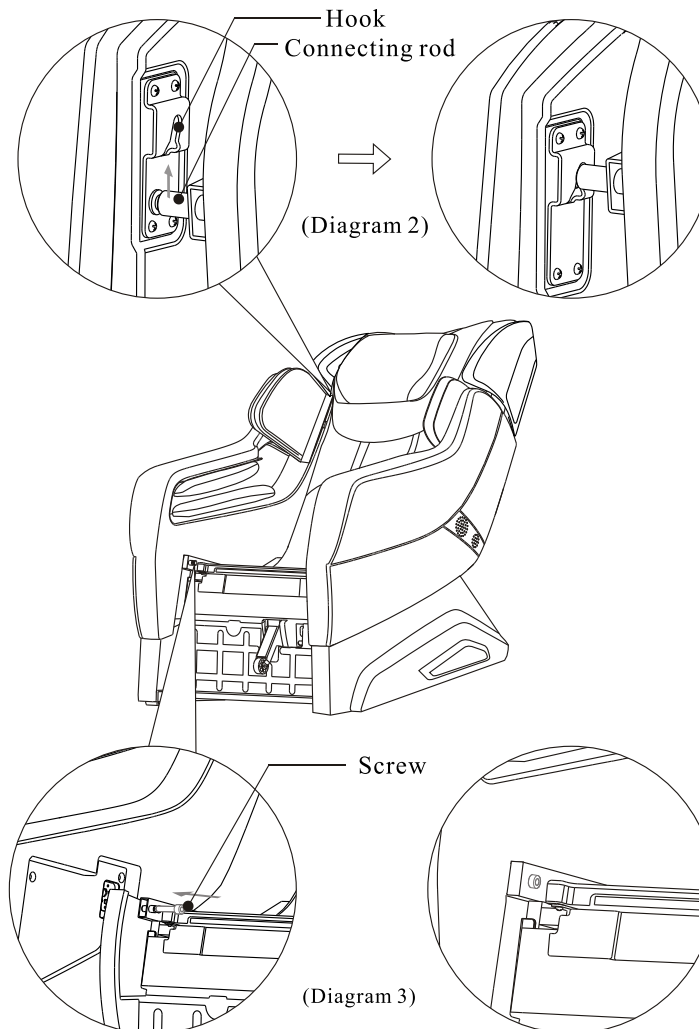

## ► PREPARATION BEFORE USING


### Installation Method

Step 1: Get out each part of the chair

- Unpack the cartons, carefully get out sewing units, accessories box and so on.
- Check all accessories with packing list.

Step 2: Installation of armrest

- Connect the power code and remote control to the relevant socket on the power box, turn on the power and make the backrest to the original position. (follow the sequence in the diagram below)
- Match the air valve, air-pipe plugs on armrest to the side ports on lower part of the armrest frame. (diagram 1)

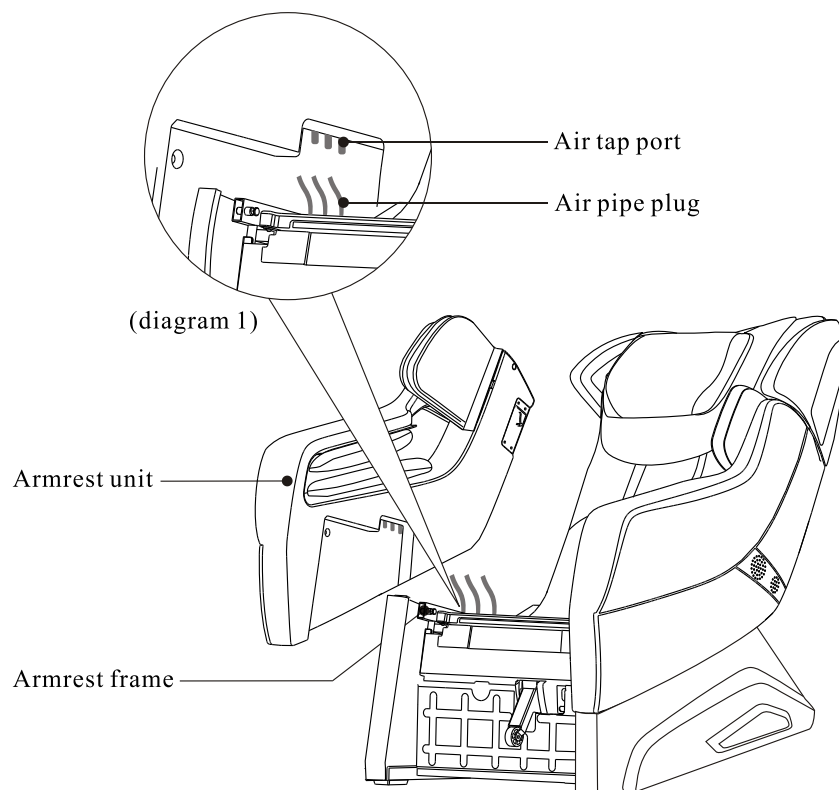


## ► NAME AND FUNCTION OF EACH COMPONENT


**PART:** Under AIRBAG PART function, press right button first, then press up or down button to choose full, arm&shoulder, back&waist, leg&foot, or buttock.

**FORCE:** Under AIRBAG FORCE function, press right button first, then press up or down button to choose airbag force from 1 to 5 grades and stop the airbag.

► **Heat:** Click "MENU", choosing "HEAT", then click right to choose on or off.

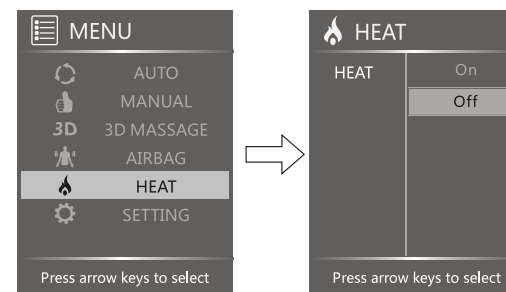

► **SETTING:** Press MENU button, choose SETTING, then press up or down button to set time, language or bluetooth.

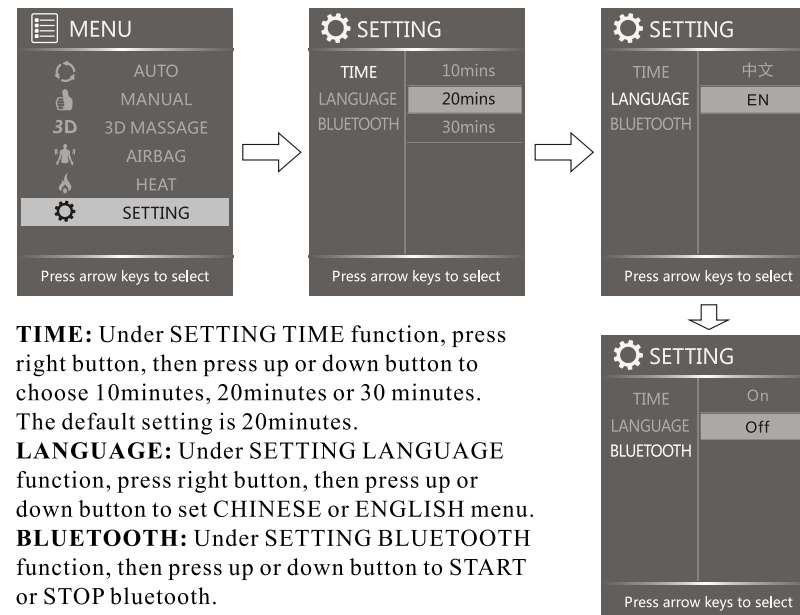

**TIME:** Under SETTING TIME function, press right button, then press up or down button to choose 10minutes, 20minutes or 30 minutes. The default setting is 20minutes.

**LANGUAGE:** Under SETTING LANGUAGE function, press right button, then press up or down button to set CHINESE or ENGLISH menu.

**BLUETOOTH:** Under SETTING BLUETOOTH function, then press up or down button to START or STOP bluetooth.


**MEMO:** All the above settings are remembered even when the machine is cut off suddenly.




## NAME AND FUNCTION OF EACH COMPONENT

Note:

- Under auto mode or manual mode, it will go back to home page as below drawing 10 seconds after without any operating.

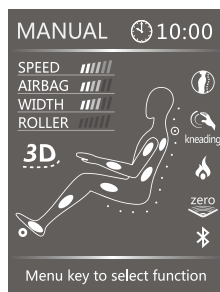

Home page in manual mode

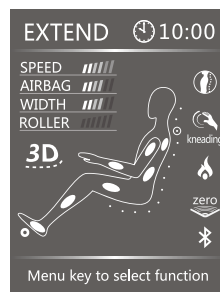

Home page in auto mode

### ⚠ Caution

- ▶ During working, click “MENU” twice, it will go back to home page immediately.
- ▶ In the home page, under manual mode or auto mode, click the left and right button to adjust air pressure intensity; under manual mode, click up and down, to adjust massage wheels speed.
- ▶ Click “AUTO” to choose auto program directly.

### ⚠ Caution

- ▶ Please do not disassemble the product to avoid any defective problems.
- ▶ Please do not splash water and other liquid into the controller to avoid malfunction of keys and error code.
- ▶ Please do not put heavy things on the controller or power cable.

## FUNCTIONS AND FEATURES

### 1. Audio Features

Product is designed with Bluetooth audio connection function.

### 2. Shoulder location

Shoulder location function gives the most effective massage for people of different height.

### 3. Zero Gravity

This chair reflects zero gravity principle from the space capsule seat of NASA, the theoretical basis of which is that the legs are higher than heart, it can imitate weightless status of aerospace, thus fully releasing the weight of body, then body will stay in most comfortable status without any stress.

### 4. 3D massage mechanism

The innovative 3D massage mechanism and traveling structure perfectly match ergonomic human engineering.

### 5. Air pressure massage

Air pressure massage for back, waist, shoulder, arm, buttock, calf, and feet.

### 6. Roller massage for sole


Double row rollers massage on foot conduct an all-around shiatsu massage onto the sole.

### 7. Diversified massage modes

Diversified massage modes: knocking, kneading/knocking, kneading, shiatsu, which act as a professional massagist.

### 8. Swing function

Product features front and rear swing function that takes you into the fantasy space swinging state and a thorough relaxation.

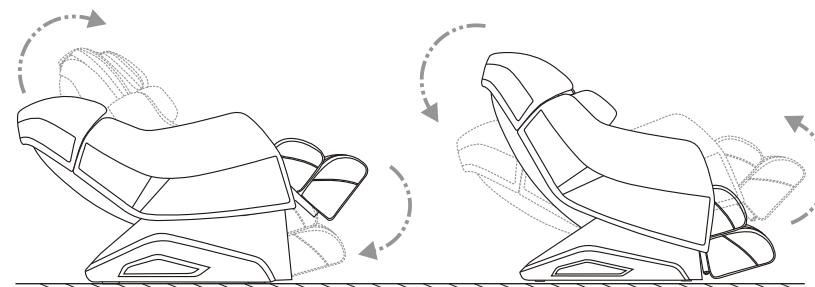

### 9. L shape super long massage track